



GESKE

GERMAN BEAUTY TECH

EDUCATION MANUAL

2023

# GESKE GERMAN BEAUTY TECH

**GESKE**  
GERMAN BEAUTY TECH

## Reimagine your skin

---

Until now, most people have taken care of their skin with cosmetic products such as creams, serums and masks. It is generally known that the best and most visible results, especially regarding pore cleansing, anti-aging treatments and lifting, are often only achieved in cosmetic salons, dermatological practices or beauty clinics. But how many people can afford that?

GESKE takes skin care to a new level.

Now the best skin care treatments like HydraFacial, MicroNeedling, LED Light Treatment and Microdermabrasion can be done at home with a handy beauty device. The price of the device is many more times cheaper than a treatment in a beauty salon.

The GESKE devices, developed together with top German dermatologists, take the skin care routine to a new level. Using beauty devices together with skin care products improves the absorption of the cosmetic active ingredients and uses advanced technologies to contain and reverse the effects of aging processes and environmental influences.





# USPs

GESKE was born after four years of in-depth research into beauty technology with the possibility to now make Artificial Intelligence (AI) available. Our range of 250+ beauty devices (and more in constant development) are a game changer.

---

### **Developed with Dermatologists**

GESKE Beauty devices have been developed in close collaboration with leading dermatologists in Germany. These state-of-the-art devices integrate advanced technologies of thermo- and cryo treatment, LED light treatment, ion treatment, electro treatment, as well as sonic pulsations to achieve visible results.

### **Clinically tested.**

To guarantee the highest quality, all GESKE Beauty devices have been extensively examined in particularly long-term application tests by the independent research institute Dermatest® and awarded the 5-star "excellent" seal.

### **Guaranteed to last for 15 years**

While a typical cream tube lasts only a few weeks, GESKE devices provide a remarkable 15-year guarantee, alongside local product warranties. Our commitment is to deliver top-notch quality at unbeatable prices.

### **Affordable**

Our devices are priced up to 90% less than market competitors, democratizing access to flawless skin for the majority of beauty consumers.

### **Turbo-charged by the GESKE app**

Our free AI app swiftly analyzes your skin with a unique Skin Scan, providing a personalized skincare routine and smart product recommendations in seconds. With 56,565 device-guided video sessions available in 45 languages, it enhances the GESKE consumer's skincare journey.

### **Attractive, easy to use, and available in 11 colors**

Beauty devices should be enjoyable. With 11 vibrant colors and user-friendly, ergonomic designs, we make skincare fun, desirable, and highly effective.

### **Most Awarded Beauty Tech Company**

Since 2022, more than 67 awards worldwide for design, app and innovation in the beauty sector.

---

We are the ONLY brand to offer such a wide range of affordable SmartAppGuided™ beauty devices and skincare products all over the world. GESKE is the obvious and logical choice for anyone who yearns to look their beautiful best.

# CLINICALLY TESTED

Developed with dermatologists.  
Clinically tested.

To ensure the highest quality, we had our products independently evaluated by the research institute Dermatest®.

In this assessment, all GESKE Beauty devices received the highest possible rating and were honored with the original Dermatest 5-star seal of excellence.

The 5-star seal is granted based on long-term application tests, providing the utmost protection to largely eliminate potential skin reactions for users in the future.

**GESKE**  
GERMAN BEAUTY TECH







# PORTFOLIO BEAUTY DEVICES

## CLEANSING



Facial Brush | 4 in 1



Facial Brush | 4 in 1 handle



Sonic Facial Brush | 5 in 1

BEST SELLER



Sonic Thermo Facial Brush | 6in1

BEST SELLER



Sonic Thermo Facial Brush & Face Lifter | 8in1

BEST SELLER



MicroCurrent Skin Scrubber & Blackhead Remover | 9 in 1

BEST SELLER

## GLOW



Sonic Warm & Cool Mask | 9 in 1



Warm & Cool Eye Energizer | 6 in 1

BEST SELLER

## LIFTING



MicroCurrent Face-Lift Pen | 6 in 1



MicroCurrent Face-Lifter | 6 in 1

BEST SELLER



Sonic Facial Roller | 4 in 1



Sonic Facial & Body Roller | 4 in 1



MicroNeedle Face & Body Roller | 8 in 1



MicroNeedle Face & Body Roller | 9 in 1



MicroNeedle Face Roller | 9 in 1

## OTHERS



Facial Hydration Refresher | 4 in 1



Facial Hydration Refresher | 4 in 1



Lip Volumizer & Booster | 4 in 1



Lip Volumizer & Booster | 4 in 1



# PORTFOLIO BEAUTY DEVICES HELLO KITTY

CLEANSING



Hello Kitty  
Facial Brush | 3 in 1



Hello Kitty  
Facial Brush | 3 in 1



Hello Kitty  
Facial Brush | 3 in 1



Hello Kitty  
Facial Brush | 4 in 1



Hello Kitty  
Sonic Thermo Facial Brush | 4 in 1



Hello Kitty  
Sonic Thermo Facial Brush | 5 in 1



Hello Kitty  
MicroCurrent Skin Scrubber &  
Blackhead Remover | 9 in 1

GLOW



Hello Kitty  
Sonic Warm & Cool Mask | 8 in 1



Hello Kitty  
Sonic Warm & Cool Mask | 8 in 1



Hello Kitty  
Warm & Cool Eye Energizer | 6 in 1

LIFTING



Hello Kitty  
MicroCurrent Face-Lift Pen | 6 in 1

OTHERS



Hello Kitty  
Facial Hydration Refresher | 4 in 1

# PORTFOLIO COSMETIC PRODUCTS

## CLEANSING



Cleansing Lotion



Facial Cleanser



Hydrating Cleanser



Cleansing Peel

## MASKS



Cleansing Mask



Ultra Hydrating Mask



Youth Boosting Mask

## SERUMS & GELS



Blue Light Protection Serum



Hydrating Serum



Hyaluronic Acid Serum



De-Puffing Eye Serum



Skin Smoothing Serum



Hydrating MicroCurrent Gel

## CREAMS



Hydrating Day Cream



UV Defense Day Cream



Anti-Aging Day Cream



Anti-Aging Night Cream







Recovery Night Cream



# ENTRY-LEVEL GESKE PRODUCTS UNDER 11€







Facial Brush | 4 in 1

-  Skin Scan & Personal Routine Guide
-  Super Soft & Flexible Silicone Tips
-  Ultra Hygienic Deep Cleansing
-  Flexible Eye-Nose Area Targeted Design



With Handle





Facial Brush | 4 in 1

-  Skin Scan & Personal Routine Guide
-  Super Soft & Flexible Silicone Tips
-  Ultra Hygienic Deep Cleansing
-  Flexible Eye-Nose Area Targeted Design



Small Size





Lip Volumizer & Booster | 4 in 1

-  Skin Scan & Personal Routine Guide
-  Natural Lip Color Revitalizing
-  Targeted Lip Volumizer
-  Ultra-Volume Technology



Big Size




Lip Volumizer & Booster | 4 in 1

-  Skin Scan & Personal Routine Guide
-  Natural Lip Color Revitalizing
-  Targeted Lip Volumizer
-  Ultra-Volume Technology

# ENTRY-LEVEL HELLO KITTY PRODUCTS UNDER 15€






### Facial Brush | 3 in 1

-  Skin Scan & Personal Routine Guide
-  Super Soft & Flexible Silicone Tips
-  Ultra Hygienic Deep Cleansing

Available Colors: ● ●




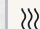

### Facial Brush | 3 in 1

-  Skin Scan & Personal Routine Guide
-  Super Soft & Flexible Silicone Tips
-  Ultra Hygienic Deep Cleansing

Available Colors: ● ●



### Facial Brush | 3 in 1

-  Skin Scan & Personal Routine Guide
-  Super Soft & Flexible Silicone Tips
-  Ultra Hygienic Deep Cleansing

Available Color: ●



**GESKE**

GERMAN BEAUTY TECH

PRODUCT TECHNOLOGIES

# SMARTSONIC PULSATION

The mild pulsations travel across the surface layers of the skin. They deliver a gentle, low-frequency massage that aids in relaxing tension points in facial and neck muscles, as well as enhancing blood circulation. In combination with the deep-warming technology the SmartSonic Pulsation Technology temporarily opens up the pores so that the active ingredients in your skin care products are able to effectively enter the stratum corneum (the outermost layer of the skin) without causing the pores to clog.

## BENEFITS

- **Ultra Hygienic Deep Cleansing**

SmartSonic Pulsations help to cleanse the skin by removing dirt, oil, and impurities more effectively than manual cleansing.

- **Exfoliation**

The SmartSonic Pulsation Technology aids in the gentle removal of dead skin cells, promoting a smoother complexion.

- **Pore Refinement**

Regular use may contribute to the appearance of smaller pores, giving the skin a smoother and more refined texture.

- **Improved Blood Circulation**

SmartSonic Pulsations can enhance blood flow, contributing to a healthier skin appearance.

- **Relaxation and Massage**

SmartSonic Pulsations offer a massaging effect, promoting relaxation and potentially reducing muscle tension in the face.

- **Stimulation of Collagen**

SmartSonic Pulsations can stimulate collagen production, which is crucial for maintaining skin elasticity.

- **Suitable for All Skin Types**

GESKE Beauty Devices allow users to adjust the intensity of SmartSonic Pulsations according to their skin sensitivity and needs.





# SMARTSONIC PULSATION DEVICES



Sonic Facial Brush | 5 in 1

- Skin Scan & Personal Routine Guide
- SmartSonic Pulsation Technology (7,500 RPM)
- Flexible Eye-Nose Area Targeted Design
- Ultra Hygienic Deep Cleansing
- Youthful Glow Massage



Sonic Thermo Facial Brush | 6 in 1

- Skin Scan & Personal Routine Guide
- SmartSonic Pulsation Technology (7,500 RPM)
- Flexible Eye-Nose Area Targeted Design
- Ultra Hygienic Deep Cleansing
- Youthful Glow Massage
- Pore-Opening Deep Warming Technology



Sonic Thermo Facial Brush & Face Lifter | 8 in 1

- Skin Scan & Personal Routine Guide
- SmartSonic Pulsation Technology (7,500 RPM)
- Flexible Eye-Nose Area Targeted Design
- Ultra Hygienic Deep Cleansing
- Youthful Glow Massage
- Pore-Opening Deep Warming Technology
- MicroCurrent Face-Lift Technology
- Ultra-Gentle Contouring & Firming Spheres


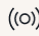


## GOOD TO KNOW

- Pulsation** With up to 18 intensities and 7,500 RPM
- Deep Warming** Steady heating temperature of the surface up to 41-43°C for pore-opening effects.

# HELLO KITTY DEVICES WITH SMARTSONIC PULSATION




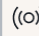


Sonic Facial Brush | 4 in 1

-  Skin Scan & Personal Routine Guide
-  Smart Sonic Pulsation Technology (7,500 RPM)
-  Ultra Hygienic Deep Cleansing
-  Youthful Glow Massage

Available Colors: ● ●




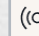



Sonic Facial Brush | 4 in 1

-  Skin Scan & Personal Routine Guide
-  Smart Sonic Pulsation Technology (7,500 RPM)
-  Ultra Hygienic Deep Cleansing
-  Youthful Glow Massage

Available Color: ●



Sonic Thermo Facial Brush | 5 in 1

-  Skin Scan & Personal Routine Guide
-  Smart Sonic Pulsation Technology (7,500 RPM)
-  Ultra Hygienic Deep Cleansing
-  Youthful Glow Massage
-  Pore-Opening Deep Warming Technology


Available Colors: ● ●

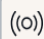


# SMARTSONIC PULSATION




Sonic Facial Roller | 4 in 1

 Skin Scan & Personal Routine Guide


 SmartSonic Pulsation Technology  
(14,000 RPM)

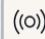
 Smoothing & Tightening Skin Massage

 3D Contouring & Firming Waves




Sonic Facial & Body Roller | 4 in 1

 Skin Scan & Personal Routine Guide

 SmartSonic Pulsation Technology  
(14,000 RPM)

 Smoothing & Tightening Skin Massage

 3D Contouring & Firming Waves

# ION TECHNOLOGY

Ion technology relies on the electrostatic interaction between ions and skin particles. Positively charged ions (cations) can attract dirt and impurities, while negatively charged ions (anions) tend to bind impurities and soothe the skin. Through targeted application, positive ions can open pores and promote the absorption of active ingredients, while negative ions relax the skin and aid in detoxification.

## BENEFITS

### ION +

#### Hygienic Skin Exfoliation

Positive ions are effective in attracting and neutralizing negatively charged impurities on the skin's surface, such as dust and pollutants. This action assists in a thorough cleansing, unclogging pores, and removing debris.

#### Detoxification

Positive ions assist in drawing out toxins from the skin, supporting the natural detoxification process and promoting overall skin health.

#### Anti-Aging Effects

Positive ions contribute to skin tightening and firming, reducing the appearance of fine lines and wrinkles. This anti-aging benefit promotes a more youthful and lifted skin appearance.

### ION -

#### Improved Hydration

Negative ions enhance the skin's ability to retain moisture, promoting improved hydration. This is beneficial for individuals with dry or dehydrated skin, leading to a softer and more supple complexion.

#### Skin Calming

Negative ions have a calming effect on the skin, reducing inflammation and soothing irritated or sensitive skin conditions. This makes them suitable for individuals with redness or skin sensitivity.

#### Stress Reduction

Negative ions are associated with stress reduction and relaxation. Incorporating negative ion technology in skincare devices can provide a soothing experience, promoting overall skin wellness through stress alleviation.





# ION TECHNOLOGY DEVICES



## MicroCurrent Skin Scrubber & Blackhead Remover | 9 in 1

### CAUTION

#### Usage

Always use with GESKE Hydrating MicroCurrent Gel or a GESKE face serum.



Skin Scan & Personal Routine Guide



Pore-Opening Deep Warming Technology



SmartSonic Pulsation Technology



ION-Deep Nurturing Technology



Red Light Active Regeneration Technology



ION+ Ultra Hygienic Deep Cleansing



MicroCurrent Face-Lift Technology



Dry and Dead Cell Removing Technology



Youthful Glow Massage

Available Colors:

Hello Kitty:

# MICRO CURRENT

Microcurrent devices activate muscle activity, thereby enhancing facial muscle tone and elevating the levels of a compound known as adenosine triphosphate (ATP), which is utilized by your cells as an energy source. ATP plays a critical role in cellular functions, including the creation of collagen. To witness the full, long-term results, it is essential to undergo a series of treatments over a period of time (ideally 30-60 days). Microcurrent has the potential to increase ATP concentration in cells by up to 400%, thus stimulating the production of collagen and elastin, effectively combating the visible signs of aging.

## BENEFITS

- **Skin Tightening and Firming**

Microcurrent treatments stimulate facial muscles, prompting them to contract and tighten. This leads to improved muscle tone, reducing sagging and imparting a firmer, lifted appearance to the skin.

- **Reduction of Wrinkles**

The stimulating effects of Microcurrent can help smooth and reduce the appearance of fine lines and wrinkles by encouraging the muscles to remain toned.

- **Collagen Production**

Microcurrent can promote the production of collagen, a protein crucial for maintaining skin elasticity. By boosting collagen levels, the treatment may support the skin's ability to bounce back and resist the formation of wrinkles and sagging.

- **Improved Blood Circulation**

The electrical currents stimulate blood circulation, enhancing the flow of oxygen and nutrients to skin cells. This improved circulation contributes to a healthier complexion, promoting cell renewal and overall skin vitality.

- **Lymphatic Drainage**

Microcurrent aids in lymphatic drainage, helping the body eliminate excess fluids and toxins. This reduces facial puffiness and swelling, resulting in a more sculpted and defined facial contour.

- **Immediate Skin Radiance**

Microcurrent sessions often provide immediate benefits, including enhanced skin radiance. The treatment boosts cellular activity and circulation, giving the skin a refreshed and revitalized appearance, making it a popular choice for pre-event beauty treatments.





# MICRO CURRENT DEVICES



MicroCurrent  
Face-Lifter | 6 in 1

- Skin Scan & Personal Routine Guide
- SmartSonic Pulsation Technology (7.500 RPM)
- MicroCurrent Face-Lift Technology
- Anti-Aging Massage Technology
- Low-Voltage Current
- Ultra-Gentle Contouring & Firming Spheres



MicroCurrent  
Face-Lift Pen | 6 in 1

- Skin Scan & Personal Routine Guide
- SmartSonic Pulsation Technology (14.000 RPM)
- MicroCurrent Face-Lift Technology
- Anti-Aging Massage Technology
- Low-Voltage Current
- Ultra-Gentle Contouring & Firming Spheres



MicroCurrent  
Face-Lift Pen | 6 in 1

- Skin Scan & Personal Routine Guide
- SmartSonic Pulsation Technology (14.000 RPM)
- MicroCurrent Face-Lift Technology
- Anti-Aging Massage Technology
- Low-Voltage Current
- Ultra-Gentle Contouring & Firming Spheres

**CAUTION**

**Usage** Always use with GESKE Hydrating MicroCurrent Gel or a GESKE face serum.

Available Colors:

Hello Kitty:

# MICRONEEDLING

Microneedling is an advanced skin treatment that involves the use of a specialized device with fine needles. These needles are gently introduced into the skin, creating tiny micro-injuries. The Microneedle device penetrates the top layer of the skin, known as the epidermis. The size and depth of these injuries can vary depending on the treatment goals. The tiny injuries signal the skin that it has been damaged, triggering a natural healing response. As a result, the body initiates repair mechanisms stimulating the skin's natural healing response. Thereby increasing the production of collagen and elastin, which are responsible for skin tightening and firmness. Through the heightened production of these substances, the skin can appear firmer and more youthful.

## BENEFITS

- **Skin Rejuvenation**

The stimulation of collagen and other growth factors can contribute to overall skin rejuvenation, providing a more youthful appearance and a fresher complexion.

- **Collagen Stimulation**

The micro-injuries created by the needles stimulate the skin's natural healing process, promoting the production of collagen and elastin. This can result in firmer and more elastic skin.

- **Improved Skin Texture**

Microneedling can enhance skin texture by reducing unevenness, fine lines, and pore size, leading to smoother and more even skin.

- **Reduced Hyperpigmentation**

Microneedling may help reduce dark spots and hyperpigmentation by promoting skin cell renewal.

- **Enhanced Absorption of Skincare Products**

The microchannels created by microneedling can improve the absorption of skincare products. After the treatment, topical products can penetrate the skin more effectively, enhancing their efficacy.

- **Minimal Downtime**

Compared to deeper microneedling treatments or other skin rejuvenation procedures, a treatment with a needle length of 0,3 mm typically has minimal downtime. Skin redness usually subsides quickly, and most people can return to their regular activities promptly.




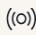



# MICRONEEDLING DEVICES



MicroNeedle  
Face & Body Roller | 8 in 1

 Skin Scan & Personal Routine Guide


 SmartSonic Pulsation Technology (14.000 RPM)


 Red Light Active Regeneration Technology

 Energizing Cooling Technology

 Impurities Prevention Technology


 DeepDerma Skin Renewal & Firming Technology

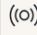
 Cell Stimulation System


 Fine Lines Refinement Technology



MicroNeedle  
Face & Body Roller | 9 in 1


 Skin Scan & Personal Routine Guide


 SmartSonic Pulsation Technology (14.000 RPM)


 Red Light Active Regeneration Technology


 Energizing Cooling Technology

 Impurities Prevention Technology

 DeepDerma Skin Renewal & Firming Technology

 Cell Stimulation System

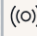
 Fine Lines Refinement Technology

 Detox Rose Quartz Spa Session




MicroNeedle  
Face Roller | 9 in 1


 Skin Scan & Personal Routine Guide


 SmartSonic Pulsation Technology (14.000 RPM)


 Red Light Active Regeneration Technology


 Energizing Cooling Technology

 Impurities Prevention Technology

 DeepDerma Skin Renewal & Firming Technology

 Cell Stimulation System

 Fine Lines Refinement Technology

 Detox Rose Quartz Spa Session

## GOOD TO KNOW

**Needles** GESKE MicroNeedle Rollers have 540 / 1.680 stainless steel needles with a length of 0,3 mm

Available Colors: 

# CRYO TREATMENT

Cryo treatment in skincare involves a non-invasive facial treatment where the skin is exposed to cold temperatures. The profound cold has both healing and revitalizing effects, as it leads to the constriction and tightening of blood vessels. Consequently, when the skin returns to its usual temperature, the blood vessels rapidly dilate, promoting improved circulation. Cryo-facials may also stimulate the production of collagen, resulting in enhanced skin elasticity.

## BENEFITS

### De-Puffing

Cooling tools help reduce facial swelling, especially under the eyes, contributing to a fresher and more alert appearance.

### Anti-Inflammatory

The cooling effect can alleviate skin inflammation and reduce redness, particularly beneficial for sensitive or irritated skin.

### Pore Refinement

By tightening the skin, cooling tools may temporarily refine pores, leading to a smoother skin texture.

### Blood Circulation

Cooling promotes enhanced blood circulation, supporting better nutrient delivery to the skin and contributing to a healthier complexion.

### Revitalization

Cooling beauty tools offer a refreshing experience that can revitalize the skin, especially in the morning or after a long day.

### Stress Reduction

The cooling massage can help alleviate tension and provide a relaxing experience, benefiting not only the skin but also overall well-being.

### Pain Relief

Cryo treatment can alleviate pain associated with certain skin conditions, injuries, or muscle soreness.





# CRYO TREATMENT DEVICES



Warm & Cool Eye Energizer | 6 in 1

Warm & Cool Eye Energizer | 6 in 1

## GOOD TO KNOW

**Deep Cooling** Steady cooling of the surface with up to  $12,5 \pm 2,5^{\circ}\text{C}$

**Deep Warming** Steady heating of the surface with up to  $39 \pm 2^{\circ}\text{C}$

Skin Scan & Personal Routine Guide

SmartSonic Pulsation Technology (11.000 RPM)

Cryo Deep Cooling Revitalizing Technology

High Efficiency Depuffing System

Youthful Glow Massage

Relaxing Eye-Warming Technology

Skin Scan & Personal Routine Guide

SmartSonic Pulsation Technology (11.000 RPM)

Cryo Deep Cooling Revitalizing Technology

High Efficiency Depuffing System

Youthful Glow Massage

Relaxing Eye-Warming Technology

Available Colors:

Hello Kitty:

# THERMO TREATMENT

Thermo treatment refers to the use of mild heat applied to the skin. This treatment is well-liked in beauty spas due to its ability to promote a feeling of overall well-being and relaxation, while simultaneously enhancing blood circulation in the skin. This improved blood flow enables oxygen and nutrients to more effectively travel through the blood vessels, a crucial process for the regeneration of damaged tissue, the activation of fibroblasts, and the facilitation of lymphatic drainage.

## BENEFITS

- **Muscle Relaxation**

The warmth can help relax facial muscles, particularly beneficial for tension and stress relief. This can contribute to a more relaxed appearance and a sense of overall well-being.

- **Improved Blood Circulation**

The warmth promotes better blood circulation in the skin, enhancing oxygen and nutrient delivery. This can contribute to a more radiant complexion and overall healthier skin.

- **Increased Absorption of Skincare Products**

The warmth opens up pores, improving the absorption of skincare products. This allows for more efficient uptake of active ingredients, enhancing the effectiveness of skincare products.

- **Stimulation of Collagen Production**

The warmth may stimulate collagen production, leading to firmer and more elastic skin. This can help to reduce the appearance of fine lines and wrinkles.

- **Skin Detoxification**

By expanding blood vessels, warmth can assist in removing toxins from the skin, contributing to improved skin health.

- **Relaxation and Stress Reduction**

The warm treatment can aid in stress reduction and provide a relaxing experience. Stress reduction is crucial for skin health, as excessive stress can exacerbate skin issues.





# LED LIGHT

NASA originally developed LED light treatment for plant growth experiments during shuttle missions. It has since shown promise in wound healing and tissue growth. LED treatment is noninvasive and works by penetrating different skin layers, using various color wavelengths to trigger natural reactions. Consistent sessions are key for noticeable skin improvements.

LED light treatment is based on the emission of light that penetrates the skin. Depending on the wavelength, this light can trigger different biological reactions in the cells. It affects the activity of cells, tissues and enzymes in a positive way.

## BENEFITS

### Red LED light

Wavelength 620-630nm

stimulates the production of collagen in the skin, thus creating a tightening and lifting effect. It also may reduce inflammation and stimulates the production of a protein responsible for younger-looking skin.

### White LED light

Wavelength = 622+470+526nm

may activate natural renewal of the skin, thus tightening it. tightens the skin and White LED light also reduces puffiness as it reaches deeper into the skin than any other LED light color.

### Blue LED light

Wavelength 465-475nm

destroys acne-causing bacteria and may help prevent breakouts.

### Green LED light

Wavelength 526nm

may help with dark circles, while also softening age spots and evening skin tone.

### Yellow LED light

Wavelength 580-590nm

may alleviate rosacea, combat swelling, soothe sunburn and improve overall skin tone.

### Cyan LED light

Wavelength 490-520nm

reduces the size of swollen capillaries, promote healing, reduces pain. Furthermore, it have an anti-inflammatory effect which calms stressed-out skin.

### Orange LED light

Wavelength 600-610nm

may revitalize the skin for a natural and healthy glow and could ease sun damage.

### Purple LED light

Wavelength 400-435nm

is a combination of blue and red LED light. This color may soothe pacify inflammation, boost radiance and support the skin to eliminate toxins.



# THERMO TREATMENT DEVICES WITH LED LIGHT



Sonic Warm & Cool  
Mask | 8 in 1

-  Pore-Opening Deep Warming Technology
-  Skin Scan & Personal Routine Guide
-  Youthful Glow Massage Technology
-  High Efficiency Depuffing System
-  Cryo Deep Cooling Revitalizing Technology
-  Full-Spectrum LED Light Technology
-  SmartSonic Pulsation Technology (8,500 RPM)
-  Full Facial Hyper- Infusion Technology

Available Colors: ● ●



Sonic Warm & Cool  
Mask | 8 in 1

-  Pore-Opening Deep Warming Technology
-  Skin Scan & Personal Routine Guide
-  Youthful Glow Massage Technology
-  High Efficiency Depuffing System
-  Cryo Deep Cooling Revitalizing Technology
-  Full-Spectrum LED Light Technology
-  SmartSonic Pulsation Technology (8,500 RPM)
-  Full Facial Hyper- Infusion Technology

Available Color: ●



Sonic Warm & Cool  
Mask | 9 in 1

-  Pore-Opening Deep Warming Technology
-  Skin Scan & Personal Routine Guide
-  Youthful Glow Massage Technology
-  High Efficiency Depuffing System
-  Cryo Deep Cooling Revitalizing Technology
-  Full-Spectrum LED Light Technology
-  SmartSonic Pulsation Technology (11,000 RPM)
-  Full Facial Hyper- Infusion Technology
-  Eye-Nose Area Targeted Design

## HOW TO USE



Take off the magnetic outer ring and place the GESKE tissue on the device's metal surface.

Afterwards apply the cleaning mask on the tissue and start the treatment.

Available Colors: ● ● ● ● ● ● ● ● ●



# HYDRATION MIST

Experience a radiant glow with our ultra-small water molecules, which efficiently penetrate your skin's surface, ensuring effective hydration and a luminous complexion. The moisturizing cooling mist technology works by splitting water droplets into ultra-small molecules in the form of mist that helps the effective components penetrate past the surface layer of the skin. Thanks to the cooling mist technology the skin will be perfectly moisturized and rejuvenated. It does not affect makeup and can be used during long travels, hot summers or just as a daily refreshment boost for the skin.

## BENEFITS

- **Hydration Boost**

Water is a fundamental element for skin health. A water spray mist can help moisturize the skin, protecting it from dehydration, particularly useful in dry or air-conditioned environments.

- **Refreshment:**

A water spray mist provides an instant feeling of refreshment, especially on warm days. This can help soothe the skin and alleviate discomfort from heat or sun exposure.

- **Skin Tone and Radiance**

Regular use can assist in maintaining a balanced skin tone, revitalizing dull-looking skin, and enhancing its natural radiance, contributing to a more vibrant complexion.

- **Makeup Setting**

A light water mist can be used to set makeup, giving it a more natural appearance. It can help fix the makeup and keep the face looking fresh.



# HYDRATION MIST DEVICES




Facial Hydration Refresher | 4 in 1


Facial Hydration Refresher | 4 in 1


Facial Hydration Refresher | 4 in 1


 Skin Scan & Personal Routine Guide


 Skin Scan & Personal Routine Guide


 Skin Scan & Personal Routine Guide


 Moisturizing Cooling Mist


 Moisturizing Cooling Mist


 Moisturizing Cooling Mist


 Energizing Hydra Refreshing Technology


 Energizing Hydra Refreshing Technology

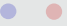
 Energizing Hydra Refreshing Technology

 3D Contouring & Firming Waves

 3D Contouring & Firming Waves

 3D Contouring & Firming Waves

Available Colors: 

Hello Kitty: 



# COSMETICS

GESKE cosmetic products feature skin-friendly ingredients and are suitable for daily use. Most of our formulas are vegan and cruelty-free. Moreover, all of the products were developed with leading German dermatologists. The cosmetic products can be paired with suitable GESKE devices for optimum results. We offer a variety of cleansers, serums, creams, and facial masks to provide the right care for every skin type.

## KEY INGREDIENTS & BENEFITS

- **Hyaluronic Acid**

Can store a lot of water, which is why it is so effective against sagging and dry skin. It is naturally produced in the body, but the older we get, the less hyaluronic acid we produce. It supports the skin's moisture supply and ensures, among other things, that wrinkles caused by dryness are reduced and the skin gains elasticity.

- **Collagen**

Collagen is a Protein produced naturally in the body, which, among other things, helps to strengthen the skin. In creams it works externally to give the skin the moisturize it needs to look fresh and radiant.

- **Vitamin B5**

Completes the cleansing effect by preventing pimples and blackheads and improving skin elasticity.

- **Vitamin E**

Promotes the innate protective barrier of the skin, increases its resistance, and contributes to a clean and fresh complexion.

- **Glycerin**

Softens the skin, plumps it up and helps retain moisture, making it a boon for dry skin.

- **Aloe Vera**

Cools and soothes the skin, moisturizes, inhibits inflammation and helps the skin to regenerate.

- **Almond Oil**

Its many unsaturated fatty acids penetrate deep into the skin and moisturize intensively. Almond Oil also contains palmitic acid, which is a natural protective component of the skin



# CLEANSERS

**GESKE**  
GERMAN BEAUTY TECH



Cleansing Lotion

For normal skin

**Vitamin E** – antioxidant, soothing  
**Almond Oil** – nourishing, moisturizing  
**Hamamelis (with Hazel)** – antioxidant, relieves itching

## KEY BENEFITS

- Intensive care and moisturization
- Non-foaming formula



Hydrating Cleanser

For normal to combination skin

**Cucumber Extract** – refreshing, anti-inflammatory  
**Provitamin B5 (Panthenol)** – moisturizing, anti-inflammatory  
**Elderflower Extract** – skin soothing

## KEY BENEFITS

- Gently removes water-based makeup
- Foaming formula



Facial Cleanser

For normal to oily skin

**Glycerin** – moisturizing  
**Provitamin B5 (Panthenol)** – moisturizing, anti-inflammatory  
**Aloe Vera** – moisturizing, refreshing, anti-inflammatory

## KEY BENEFITS

- Prepares skin for serums and moisturizers
- Foaming formula



Cleansing Peel

For normal to oily skin

**Bamboo Peeling Grains** – gently cleansing  
**Provitamin B5 (Panthenol)** – moisturizing, anti-inflammatory  
**Aloe Vera** – moisturizing, Refreshing, anti-inflammatory

## KEY BENEFITS

- Purifying facial scrub for a clear complexion
- Prepares skin for serums and moisturizers



# FACE & EYE SERUMS



Blue Light Protection Serum

For all skin types

**Malachite Extract** – protects skin from environmental influences, antioxidant  
**Laminaria Digitata** – protects the skin from external influences  
**Ginseng** – antioxidant

## KEY BENEFITS

- Protects skin from blue light damage



Hyaluronic Acid Serum

For all skin types

**Hyaluronic Acid** – smoothing, toning  
**Provitamin B5 (Panthenol)** – moisturizing, anti-inflammatory  
**Betaine** – nourishing, moisturizing

## KEY BENEFITS

- Deep moisturization
- Triple structure hyaluronic acid for maximum impact



Hydrating Serum

For normal to oily skin

**Hyaluronan** – smoothing, toning  
**Vitamin A** – skin smoothing, nourishing  
**Vitamin C + E** – antioxidant, nourishing

## KEY BENEFITS

- Radiant skin thanks to intensive care
- For a fresh glow



De-Puffing Eye Serum

For all skin types

**Hyaluronan** – smoothing, toning  
**Provitamin B5 (Panthenol)** – moisturizing, anti-inflammatory  
**Aloe Vera** – moisturizing, refreshing

## KEY BENEFITS

- Moisturizes the delicate skin around the eyes
- Removes puffiness and dark circles

# FACE MASKS



Cleansing Mask

For all skin types

**Soy Protein** – hydrating, soothing  
**Argan Oil** – antioxidant  
**Sesame Seed & Black Cumin Oil** –  
nourishing, moisturizing

## KEY BENEFITS

- Detox for the skin
- Skin refining formula



Youth Boosting Mask

For all skin types

**Hyaluronan** – smoothing, toning  
**Aloe Vera** – moisturizing, refreshing  
**Vitamins E + C** – antioxidant, nourishing

## KEY BENEFITS

- For a delicate glow



Ultra Hydrating Mask

For dehydrated skin

**Hyaluronan** – smoothing, toning  
**Shea Butter** – nourishing, moisturizing  
**Green Tea** – soothing

## KEY BENEFITS

- Glowing skin formula
- For a fresh complexion and soft skin

FOR BEST RESULTS USE WITH





# DAY CREAMS

**GESKE**  
GERMAN BEAUTY TECH



Hydrating Day Cream

For all skin types

**Almond Oil** – nourishing, moisturizing  
**Provitamin B5 (Panthenol)** – moisturizing, anti-inflammatory  
**Quince Extract** – antioxidant  
**Jojoba Oil** – calms the skin

#### KEY BENEFITS

- 24-hour moisturization and long-lasting glow
- From pale to glowing skin



Anti-Aging Day Cream

For normal to dehydrated skin

**Collagen** – skin smoothing, skin toning  
**Hyaluronan** – smoothing, toning  
**Rambutan Leaves Extract** – hydrating, smoothing

#### KEY BENEFITS

- For youthful looking skin
- SPF 20



UV Defense Day Cream

For normal to dehydrated skin

**Q10** – hydrating, smoothing  
**Collagen Complex** – skin smoothing, skin toning  
**Hyaluronan** – smoothing, toning

#### KEY BENEFITS

- Has an antioxidant effect and helps to improve cell health
- Reduces wrinkles caused by dryness

# NIGHT CREAMS



## Anti-Aging Night Cream

For normal to dehydrated skin

**Collagen** – skin smoothing, skin toning  
**Hyaluronan** – smoothing, toning  
**Rambutan Leaves Extract** – hydrating, smoothing

### KEY BENEFITS

- For youthful looking skin



## Recovery Night Cream

For normal to dehydrated skin

**Q10** – hydrating, smoothing  
**Collagen Complex** – skin smoothing, skin toning  
**Argan Oil** – moisturizing, anti-bacterial

### KEY BENEFITS

- For smoother and firmer skin



# GEL & BODY SERUM



Hydrating MicroCurrent Gel

For all skin types

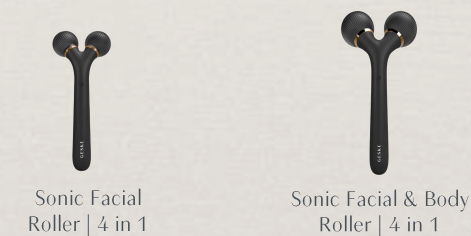
BEST RESULTS WITH



Skin Smoothing Serum

For all skin types

BEST RESULTS WITH



**Collagen** – skin smoothing, skin toning  
**Hyaluronan** – smoothing, toning  
**Glycerin** – moisturizing

## KEY BENEFITS

- Always apply before starting your MicroCurrent Session
- Moisturizes and soothes the skin

**Collagen** – skin smoothing, skin toning  
**Caffeine** – nourishing, detoxifying  
**Pomelo Oil** – nourishing, anti-bacterial

## KEY BENEFITS

- Strengthens your connective tissue
- Firming formula





# CONTACT

## E-Mail:

[trademarketing@geske.com](mailto:trademarketing@geske.com)

[patricia.sabando@geske.com](mailto:patricia.sabando@geske.com)

## Social Media:

[TikTok](#)

[YouTube](#)

[Facebook](#)

[Instagram](#)

## Address:

Leipziger Platz 18

10117 Berlin

Germany

