
Upute za upotrebu

BEEABeauty

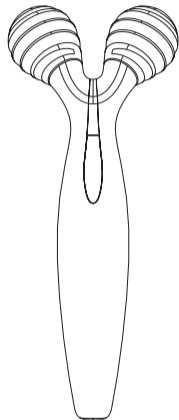TB-1682

Mikrostruja, solar panel
VALJAK ZA MASAŽU LICA

※ Da biste osigurali pravilan rad ovog uređaja, pažljivo pročitajte ove upute prije upotrebe uređaja i sačuvajte ih za buduću referencu.

UPOZORENJE

-  Ovaj uređaj nije igračka, držite ga izvan dohvata djece.
-  Molimo provjerite uređaj prije početka upotrebe. Ne koristite ga ako je oštećen.
-  Ako se uređaj koristi u suhom okruženju, rukovanje uređajem može generirati statički elektricitet. To je normalna pojava i ne utječe na upotrebu uređaja

 **OPREZ**

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
BEEABeauty

TB-1682

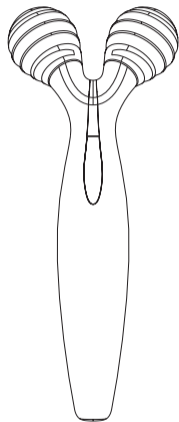Microcurrent, solar panel
FACIAL ROLLER

※ To ensure proper operation of this device, please read this instruction manual carefully before using the device and keep it for further reference.

WARNING

-  This device is not a toy, keep out of reach of children.
-  Please verify the device before you start using it. Do not use if there is any damage.
-  If used in a dry environment, handling this device may generate static electricity. This is a normal phenomenon. It does not affect the use of the device.

 **CAUTION**

 The following people should not use this device: pregnant women, while nursing or during menstruation; people with allergies; people with facial nerve disorders; people who have fractures or metal, silicone, or plastic implants; people with open wounds, chapped and swollen skin on the area of use; people with other severe skin diseases. This device is not intended for use by people (including children) with reduced physical capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning the use of the device by a person responsible for their safety. Immediately stop using the device and consult a doctor if you feel any discomfort. Keep the device out of reach of children. Regularly clean the device to extend its lifetime. Do not share the same device with others. Do not attempt to modify, disassemble, or repair the device by yourself to avoid accidents. Remove any jewelry/accessories and tie your hair when using the device, to prevent hair and jewelry from interfering with the device. Make sure the solar panel is exposed to the light source when using.

Device Specifications

TB-1682


- ① Name: Microcurrent facial roller BEEA Beauty solarni panel
- ② Model No.: TB-1682
- ③ Solar Panel Rated Voltage: 200lux light \geq 1.35V
- ④ Microcurrent: $\leq 10\mu\text{A}$

Device Description

TB-1682

- ① Massage Rollers
- ② Solar Panel
- ③ Handle
- ④ Cleaning Cloth
- ⑤ Storage Bag

Picture 1

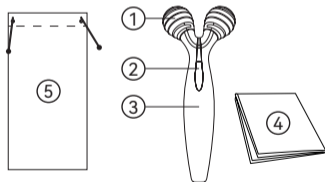Preparation Instructions:

TB-1682

※ This device does not require charging. The device features a solar panel on the handle, which uses light to generate microcurrent and output on the rollers. Make sure you do not cover the solar panel while using the device.

How to Use:

※ Clean the device with a clean damp cloth. Use the device on the skin in different angles to increase the massage intensity and efficiency. (Picture 2)

Face: Gently move the rollers from the jawline to your cheekbone, or from the corner of the mouth to the ear; move the rollers in upward movements and from the center to the sides of the face. (Picture 3)

Neck: Gently move the rollers along the neck to the jawline. (Picture 4)

Arms: Gently move the rollers along the arm. (Picture 5)

※ Do not use this device excessively. It is recommended to massage each area 15–20 times.

Picture 2

Picture 3

Picture 4

Picture 5

Cleaning & Maintenance

TB-1682

- Do not immerse the device in water. Clean the surface of the device with a clean damp cloth.
- Never use steel wool, abrasive cleaners or corrosive liquids to clean the device, may result in its malfunction or discoloration.
- Never use lubricants, detergents, polishing agents or air fresheners on this device.
- After cleaning, please use a towel to wipe the device.

Storage Do not place or store the device for a long period of time under a hot or humid environment.

Recycling When disposing of this device, do not put it in the general household waste. Please dispose of this device at a recycling point to contribute to environmental protection.

Troubleshooting

TB-1682

Problem	Possible Cause	Solution
Rollers rotation is not smooth.	Rollers have debris.	Consult your local retailer.
Other problems.	/	Consult your local retailer.

Note: If the device still does not work properly after using the above solutions, please stop using the device and contact your local retailer.

Warranty

- ※ The warranty period starts from the purchase date.
- ※ The warranty period is 1 year.
- ※ In case of failure, you can return device directly in the store where it was purchased with the purchase invoice, this warranty card, and a detailed description of the defect.

The warranty will be void under the following conditions:

- If the customer does not present the correct warranty card and purchase receipt.
- If the customer did not follow the instructions for using the product
- If the product has been opened, modified or repaired by an unauthorized person.
- If product defects were caused by force majeure such as: lightning strike, electric shock in the electrical network, natural disasters, etc.
- If the defects were caused by damage, improper use or improper transportation.
- If the failure was caused by an error in the system to which the product is connected

This guarantee does not change the legal consumer rights valid in EU in relation to those prescribed by the manufacturer.

Authorized service:

Feniks-Info d.o.o.,
Slavonska avenija 72, Zagreb,
servis@feniks-info.hr
tel 01/6395 281

Date of purchase

Stamp and signature