
Upute za upotrebu

BEEABeauty

TB-1712

Radio frekvencija, Mikrostruja, LED
ANTI AGING

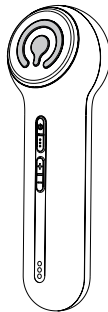⚠ UPOZORENJE

- Ovaj uređaj nije igračka, čuvajte ga izvan dohvata djece.
- Molimo provjerite uređaj prije početka uporabe. Ne koristite ga ako je oštećen.
- Ako se uređaj koristi u suhom okruženju, rukovanje uređajem može stvarati statički elektricitet. To je normalna pojava i ne utječe na uporabu uređaja.
- Molimo, držite USB kabel za punjenje i priključak za punjenje uređaja suhim.
- Ne koristite uređaj dok se puni.

⚠ OPREZ

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
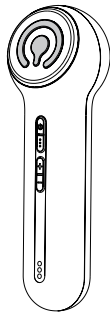
Beauty Device
ANTI AGING

WARNING

-  This device is not a toy, keep out of reach of children.
-  Please verify the device before you start using it. Do not use it if there is any damage.
-  If used in a dry environment, handling this device may generate static electricity. This is a normal phenomenon. It does not affect the use of the device.
-  Please, keep the USB charging cable and the device charging port dry.
-  Do not use the device while charging.

CAUTION

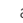-  Do not use this device in the following situations:
 -  Pregnant women or while nursing; children; people who have a pacemaker or other electronic medical equipment; people with allergic reactions; while wearing makeup.
 -  If the treatment area has skin conditions, such as open wounds, ringworm, eczema, rash, burns, sunburns, abnormal wound healing, etc.
 -  People with diabetes, lupus erythematosus, porphyria or congestive heart disease, bleeding disorder, history of immunosuppressive disease (including HIV infection and AIDS).
 -  Do not use this device in conjunction with any medical electronic equipment such as a pacemaker, internal defibrillator, etc. Do not use if you have arrhythmia or any known heart disease.
 -  People with metal implants, or piercings on the area of use.

- Do not use this device if you are receiving radiation therapy or chemotherapy treatments.
- People with diseases that can be stimulated by heat, such as herpes simplex virus.
- Do not use this device if you take any of the following medication:
 - If your skin is currently being treated with or has recently been treated with Alpha-Hydroxy Acids (AHAs), Beta-Hydroxy Acids (BHAs), topical isotretinoin and azelaic acid.
 - If you have taken any form of Accutane A or Roaccutane A within the past six months.
 - If you are taking photosensitizing agents or any other doctor- prescribed medication, please check the patient information leaflet before you start using this device. If the patient information leaflet states that it can cause photo-allergic reactions, photo-toxic reactions, or if it specifies to avoid exposure to sunlight when taking the medicine, please do not use the device.
 - If you take anticoagulant medication, including large doses of aspirin, only use the device one week after you stop the treatment.
- This device is suitable to use on the face (forehead, cheeks, chin). Never use it on the eyelids, lips, chest, genital area, throat, or on areas where you have had surgery or epilation.
 - If your skin is light sensitive or prone to allergic reactions, do not use the light therapy mode.
 - Do not use the light therapy mode if the treatment area has infections, eczema, burns, inflammation of hair follicles, open lacerations, herpes, wounds or lesions, and hematomas.
 - Do not look directly at the LED light. Use eye protection during the light therapy treatment, such as sunglasses or protective goggles (not included).

- Do not repair, disassemble or modify the device by yourself, as improper handling may damage it. Do not use this device for any other purpose.
- If you have any doubt or concern about side effects, please, consult with your doctor.
- If it is recommended to use a gel that is compatible with this device (purchase separately). If you use any other cream or gel, it may affect the normal operation of this device and RF function.
- This device has a lithium battery. Do not place or store the device for a long period in a hot or humid environment. If you do not intend to use the device for a long time, store it after it is fully charged. Do not throw the device into a fire or put it in the general household waste. Please dispose of this device in a recycling point.

Device Specifications


TB-1712

- Name: BEEA Beauty Multi-Therapy Beauty Device
- Model No.:TB-1712
- Input: 5V—/1.5A
- Power: 7.4V—/4W
- Type of Battery: Li-ion 7.4V 1000mAh

Device Description

TB-1712

- ① Electrode Poles
- ② LED Lamps
- ③ Power Button
- ④ Mode Switch Button
- ⑤ Intensity Button +
- ⑥ Intensity Button -
- ⑦ Intensity Indicator
- ⑧ RF Indicator
- ⑨ MC Indicator
- ⑩ LED Indicator
- ⑪ USB Plug
- ⑫ Contact Charging Points
- ⑬ Charging Base
- ⑭ Charging Base Contact Points
- ⑮ Charging Base Cable Port
- ⑯ Back Adhesive
- ⑰ Suspension Holes

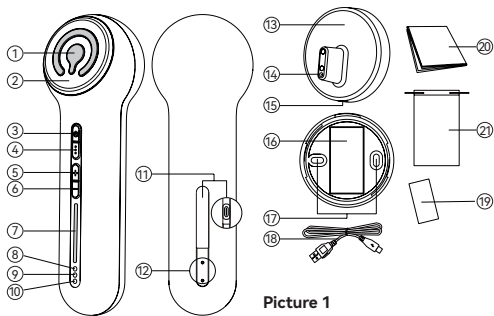
Picture 1


- ⑱ USB Cable
- ⑲ Adhesive Strip
- ⑳ Cleaning Cloth
- ㉑ Protective Bag

☹ Charging the Device

- ※ This device provides two charging systems: USB charging cable (picture 2) or wall charging base. (Picture 3). Charge the device according to the charging system. Charging the device takes approximately 3 hours. When the battery is fully charged, the device can be used consecutively for about 60 minutes. The intensity indicator blinks in a cycle when charging, and it lights up completely when the battery is full.
- ※ If the device battery is low, the intensity indicator blinks and the device automatically shuts down. The device does not work while it's charging.

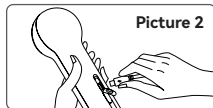
Picture 2

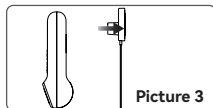Picture 3

☹ Install the Charging Base

- ※ Stick the adhesive strip to the back case, with UP sign facing up, fix it on the wall. (Picture 4)

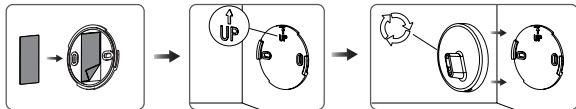Picture 4

Button Instructions

Power Button:

Press the power button. After a beeping sound, the device turns on the RF mode first intensity level. --> Press again, the device turns off with a beeping sound.

Mode Switch Button:

The device turns on the RF mode with a beeping sound. --> Press the mode switch button to change to microcurrent mode. --> Press again to change to LED mode --> cycle.

Intensity Button +: Press to increase the intensity.

Intensity Button -: Press to decrease the intensity.

Vrijeme korištenja

RF Every 60s the device emits a beeping sound indicating to change the treatment area, after 4m it automatically turns off.

RF+MC
Every 90s the device emits a beeping sound indicating to change the treatment area after 6m it automatically turns off.

RF+MC+LED Every 120s the device emits a sound indicating to change the treatment area, after 8m it automatically turns off.

How to Use

TB-1712


To ensure the correct operation of this device, please read the instructions carefully before use.

Recommended Use

- ※ Use the device on clean skin. Be sure to remove makeup, excess oil, and impurities before use.
 - ※ When you use the device, do not look directly at the LED light source to avoid discomfort.
- 1) Remove all jewelry from your body before using the device, including earrings, rings, necklaces, etc.
 - 2) Divide the treatment area into five sections: forehead, left and right cheeks, chin, and neck. (Picture 5)
 - 3) How to apply the gel: Apply enough gel or essence to cover the electrode poles surface. Place the device on your skin and move it slowly and evenly from the center of the face outwards. (Picture 6)
 - 4) Turn on the device, select the mode and intensity level you want to use. Slowly move it in a circular motion to warm up the skin.
 - 5) If the warm feeling is too low, reduce the range of the circle and speed movement; if the warm feeling is too high, increase the range of the circle and speed movement. Adjust the temperature to a comfortable level. Every 90s, the device emits a sound reminding you to change the treatment area. (Picture 7)

Slika 5

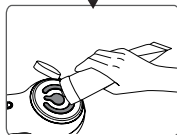
Slika 6


- 6) When changing treatment area, apply enough gel or essence to cover the electrode poles surface and repeat the above steps.
- 7) For a lifting effect, move the treatment surface from the cheek towards the ear; move upward from the eyebrows to the hairline to lift the forehead. (Picture 8)
- 8) Be sure to constantly move the device across your skin while using it. Although this device features a protection system, different people have different skin tolerance and sensitivity. Keeping the device without movement on a specific area may cause minor burns.
- 9) When using the light therapy mode, do not look at the light therapy source. (Picture 9)
- 10) After use, wipe off any remaining gel from the skin and follow your skincare routine.

Recommended Use Frequency

- It is recommended to use one mode only per session. Allow one day rest between treatments.
- It is recommended to use the device 2-3 times a week for 4-8 minutes. Use the device for a minimum of 6-8 weeks, after this time, use according to your skin needs. Use the device consistently for better results.

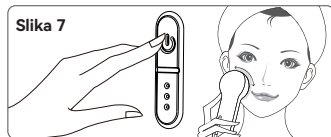Slika 8

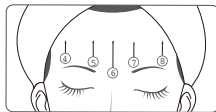Slika 9

WARNING

- 1) If you have purchased the gel, use it before the expiration date specified on the packaging. The gel is for external use only.
- 2) For first-time users, test the gel or the essence on the skin 24h prior to initial use. Apply a small amount of gel or essence on the inside of the arm to ensure it will not cause a skin reaction.
- 3) Use the device with the gel or the essence on the inside of the arm (8 min) 24h before you use the device. After use, the skin may show signs of redness or appear slightly swollen for a short time. If these reactions do not disappear within 24 hours, please consult a doctor.
- 4) Avoid using the device for too long on the cheeks or on areas you want to avoid reduce volume.
- 5) Skin sensation and results may vary from person to person and area of use.
- 6) This device has an overcurrent protection function. When the current exceeds the maximum capacity, the device will shut down automatically after three "DIDIDI" sounds.
- 7) During use, do not leave the device on one place without movement for longer than 2-3 seconds. After 2-3 seconds, make sure you start moving the device across your skin to avoid unnecessary minor burns.

Cleaning & Maintenance

TB-1712

-  After using the device, clean the electrode poles surface and fuselage with a clean slightly wet cloth. Do not submerge the device in water or any other liquid substances.
-  Do not use steel wool, abrasive cleansers, or corrosive liquids to clean the device.

Storage

This device has a lithium battery. Do not place or store the device for long periods under a hot or humid environment. If you do not intend to use the device for a long time, store it after it is fully charged.

Recycling

The built-in battery contains substances that may pollute the environment. When disposing of this device, do not put it in the general household waste. Please dispose of this device in a recycling point to contribute to environmental protection.

Troubleshooting

TB-1712

| Problem | Possible Cause | Solution |
|---|--|--|
| The device does not work properly. | Battery is low or out of power. | Recharge the device as guided. |
| Cannot feel heat or microcurrent. | Not enough gel or essence. | Apply more gel or essence. |
| | You are moving the treatment head too quickly. | Move the treatment head slowly. |
| | The device is not working. | Move the treatment head slowly and evenly. |
| Other problems | | Contact your local retailer. |
| Note: If the device still does not work properly after using the above solutions, please stop using the device and contact your local retailer. | | |

Warranty

- ※ The warranty period starts from the purchase date.
- ※ The warranty period is 1 year.
- ※ In case of failure, you can return device directly in the store where it was purchased with the purchase invoice, this warranty card, and a detailed description of the defect.

The warranty will be void under the following conditions:

- If the customer does not present the correct warranty card and purchase receipt.
- If the customer did not follow the instructions for using the product
- If the product has been opened, modified or repaired by an unauthorized person.
- If product defects were caused by force majeure such as: lightning strike, electric shock in the electrical network, natural disasters, etc.
- If the defects were caused by damage, improper use or improper transportation.
- If the failure was caused by an error in the system to which the product is connected

This guarantee does not change the legal consumer rights valid in EU in relation to those prescribed by the manufacturer.

Ovlašteni servis: Feniks-Info d.o.o.
Slavonska avenija 72, Zagreb,
servis@feniks-info.hr; tel 01/6395 281

Datum kupnje:

Pečat i potpis: