
Upute za upotrebu

BEEABeauty

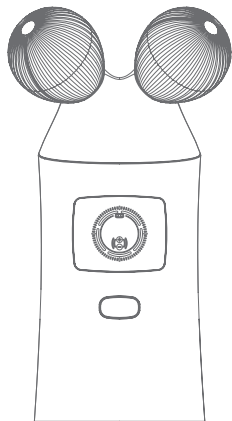TB-1767

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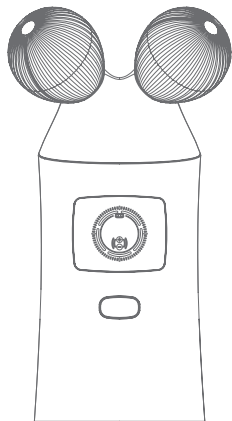
TB-1767

3-in-1
FACIAL BEAUTY DEVICE

✘ To ensure the proper operation of this device, please read this instruction manual carefully before using the device and keep it for further reference.

WARNING

-  This device is not a toy, keep out of reach of children.
-  Please verify the device before you start using it. Do not use if there is any damage.
-  Please, keep the USB cable and the device charging port dry.
-  Do not use the device while charging.

CAUTION

-  This device should not be used by:
 - Pregnant women; children; people with allergies;
 - People with skin conditions or diseases (psoriasis, burns, scalds, sunburns, etc.); people with diabetes, lupus erythematosus, porphyria or congestive heart disease, or if you have a history of immunosuppressive diseases (including HIV or AIDS);
 - People receiving radiation or chemotherapy treatments.

- ⦿ This device is not intended for use by people (including children) with reduced physical, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning the use of the device by a person responsible for their safety.
- ⦿ Do not use this device in the following areas:
 - Eyelids and throat;
 - Epilated areas; areas where surgery has been performed.
- ⦿ Do not use this device if you take any of the following medication:
 - The skin is or has been treated with alpha-hydroxy acids (AHAs), beta hydroxy acids (BHAs), local isotretinoin or azelaic acid.
 - If you have taken over the past six months oral or topic Accutane A or Roaccutane A. This treatment makes skin susceptible to lacerations, trauma, and irritation.
 - If you take photosensitizer or any other medicine, always verify its informative leaflet, if there is any reference that it may cause light allergy or phototoxicity, or if it mentions not to sunbath
do not use the light therapy head.
 - If you take anticoagulant medication, including the use of aspirin in large quantities, use the device only one week after you finish taking the treatment.
- ⦿ If your skin is sensitive to light, prone to rash or allergic reactions, do not use the light therapy.

- Do not use the light therapy head if the area you want to treat has an infection, eczema, burn, hair follicle inflammation, open wounds, scratches, herpes, lesions, or hematomas.
- Stop using the device and consult a doctor if you feel any discomfort during or after use.
- Do not use the light therapy head for a long time in the same area.
Follow this instruction manual information for precise use.
- If you are using the device for the first time, use the low-intensity mode for a better experience. Use microcurrent massage heads with a water-soluble gel for better results.
- Do not look directly at the light therapy source. It is recommended to wear sunglasses/ protective goggles (not included) to protect the eyes.
- Keep the device out of reach of children.
- Regularly clean the device to extend its lifetime.
- Do not share the same device with others.
- Do not attempt to disassemble, repair, or modify the device by yourself to avoid accidents.

Device Specifications


TB-1767

- Ⓜ Name: BEEA Beauty 3-in-1 Facial Beauty Device
- Ⓜ Model No.: TB-1767
- Ⓜ Input: 5V /0.5A
- Ⓜ Power: 3.7V /1.7W
- Ⓜ Type of Battery: Li-ion 3.7V 800mAh

Device Description

TB-1767

- ① Microcurrent Facial Massage Head
- ② Display Area
- ③ Power/Mode/Intensity Switch Button
- ④ Charging Port
- ⑤ Microcurrent Eye & Lip Massage Head
- ⑥ Light Therapy Head
- ⑦ Time Indicator
- ⑧ Battery Indicator
- ⑨ Light Therapy Mode Indicator
- ⑩ Intensity Indicator
- ⑪ Microcurrent Mode Indicator
- ⑫ USB Cable

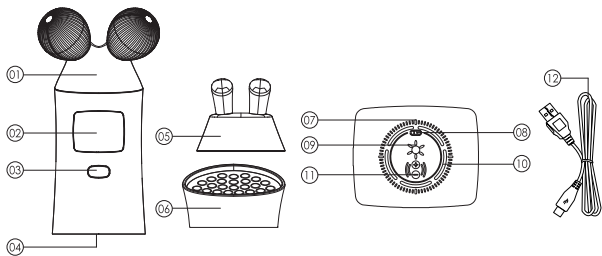PICTURE 1.

Preparing the Device for Use

TB-1767

- ※ The battery may be low before you use the device for the first time. Please charge the device with the USB cable.
- ※ The device has a low battery reminder. When the battery is low, the battery indicator will flash for 10 seconds, and the device shuts off automatically. Please charge the device.

Charging Instructions

- Make sure the device is turned off.
- Connect the USB cable to the device's USB charging port located on the bottom of the device.
- Plug the adapter into the wall socket.
- While charging, the battery indicator lights up in a cycle. The battery indicator fully lights up when the battery is charged. Charging the device takes approximately 3 hours.
- A full charge allows up to 70 minutes of use.

Caution

TB-1767


1. Do not attempt to replace the device's battery.
 2. Always use the original USB cable to charge the battery.
- ※ Use the device on clean and dry skin.
 - ※ When using the light therapy head, wear protective goggles (not included) to protect the eyes.

How to Use


TB-1767

Remove and Install the Heads:

- Remove: Hold the device's body in one hand, with the other hand rotate the head to the right and pull to remove it. (Picture 2)
- Install: Align the head with the device's body and connection buckle, gently press the head and rotate it to the left to snap in place. (Picture 3)

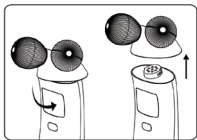
PICTURE 2.

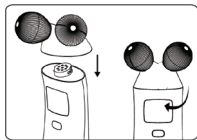PICTURE 3.

Using the Microcurrent Facial Massage Head

- ✘ This microcurrent mode has five adjustable intensity levels to fit different skincare needs.
- ✘ Use the device on the skin at different angles to increase the massage intensity and efficiency. (Picture 4)

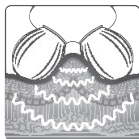PICTURE 4.

- ※ This head attachment provides a 6-minute microcurrent treatment. Every 2 minutes the device emits a “didi” sound reminding you to change the area of use.
- Install the microcurrent facial massage head.
 - Apply water-soluble gel to the areas you want to use the device.
 - Press and hold the power/mode/intensity switch button to turn on the device.
 - Press the power/mode/intensity switch button to adjust the microcurrent intensity.
 - If you are using this device for the first time, we recommended that you start at a low-intensity.
 - Cheeks: Slowly glide the massage rollers from the jawline to your cheekbone or from the corner of your mouth to the ear. We recommended to use this head for 2 minutes on each cheek. (Picture 5)
 - Forehead: Close your eyes and slowly glide the rollers from the top of the eyebrows to the hairline. It is recommended to use the device for a total of 2 minutes on the forehead. (Picture 6)

PICTURE 5.

PICTURE 6.

- Microcurrent massage helps to tone muscles, improve the appearance of wrinkles and fine lines on the chin, cheeks, and forehead. Use the device everyday for better results, or use it according to your skin needs for 6 minutes each session.
- Using the Microcurrent Eye & Lip Massage Head
 - ✘ This microcurrent mode has three adjustable intensity levels to fit different skincare needs.
 - ✘ Press the power/mode/intensity switch button to adjust the microcurrent intensity. First-time users should choose a low-intensity level.
 - ✘ This head attachment provides a 3-minute microcurrent treatment. Every minute, the device emits a “didi” sound reminding you to change the area of use.
- Install the microcurrent eye & lip massage head.
- Apply water-soluble gel to the areas you want to use the device.
- Press and hold the power/mode/intensity switch button to turn on the device.
- Press the power/mode/intensity switch button to adjust the microcurrent intensity.
- Around the Mouth: Place the massage head on the corner of your mouth and slowly glide it towards the ear. Repeat this step 3 times. It is recommended to use for 1 minute. (Picture 7)

PICTURE 7.

- Around the Eyes: Place the massage head on the corner of the eye and slowly glide or press towards the ear. Repeat this step 3 times. We recommended to use the device for 1 minute on each eye. (Picture 8)
- Microcurrent eye & lip massage head helps to reduce wrinkles and fine lines around the eyes and mouth. Use the device every day for optimal results or use it according to your skin needs for a total of 3 minutes each session.

PICTURE 8.

👁 Using the Light Therapy Head

- ✘ Before using the light therapy head, test it on your skin to make sure you are not sensitive to light. Use the light therapy head on the forearm for 3 minutes. If your skin shows signs of redness, and it lasts for more than 2 hours, your skin may be sensitive to light. In this case, do not use the light therapy head.
- ✘ When using the light therapy head, always wear goggles to avoid looking directly at the light source to protect the eyes and prevent discomfort.
- ✘ When using the light therapy head, place it in direct contact with your skin over the wanted area.

- Install the light therapy head.
- Long press the power/mode/intensity switch button to turn on the device. Press the power/mode/intensity switch button to change to light therapy mode.
- The red and yellow light modes work for 6 minutes. Every 2 minutes, the device emits a sound reminding you to change the treatment area. The blue light mode works for 3 minutes. Every minute, the device emits a sound reminding you to change the treatment area. When the treatment time ends, the device automatically turns off. (Picture 9)

PICTURE 9.

- Use the mixed-light mode for comprehensive care, including areas with acne. The mixed-light mode works for 6 minutes. Every 2 minutes, the device emits a sound reminding you to change the treatment area. The device automatically turns off after 6 minutes for precise use.
- Change the light mode according to your skin needs. Do not overuse light therapy on your skin, and always follow the recommended using time. After use, you can follow your normal skincare regimen.
- It is recommended to use the light therapy head 2-3 days a week for 6 minutes each session.

Cleaning & Maintenance

TB-1767

- ※ To ensure the best performance of the device, clean the device body and massage heads with a clean damp cloth before and after each use.
- 🚫 Make sure the device is power off before you start cleaning.
- 🚫 Do not immerse the device in water.
- 🚫 Do not use steel wool, abrasive cleansers, chemicals, or corrosive liquids to clean the device.
- 🚫 Do not use lubricants, detergents, polishing agents, or air fresheners on this device.

Storage

TB-1767

This device has a lithium battery. Do not place or store the device for long periods in a hot or humid environment. If you do not intend to use the device for a long time, store it after it is fully charged.

Recycling

TB-1767

The built-in battery contains substances that may pollute the environment. When disposing of this device, do not put it in the general household waste. Please dispose of this device at an official recycling point to contribute to environmental protection.

 Troubleshooting

Problem	Possible Cause	Solution
The device does not turn on.	Massage head was installed incorrectly.	Install the massage head correctly as guided.
	The battery is low or out of power.	Charge the device as guided.
Microcurrent sensation is weak. Massage head rotation is not smooth.	The microcurrent intensity level is low.	Follow the instructions to adjust the microcurrent intensity level.
	The skin is too dry.	Replenish your skin with water soluble-gel before you use the device.
Other problems.		Contact your local retailer.
Note: If the device still does not work properly after using the above solutions, please stop using the device and contact your local retailer.		

The warranty period starts from the purchase date.

The warranty period is 1 year.

In case of failure, you can return device directly in the store where it was purchased with the purchase invoice, this warranty card, and a detailed description of the defect.

The warranty will be void under the following conditions:

- If the customer does not present the correct warranty card and purchase receipt.
- If the customer did not follow the instructions for using the product
- If the product has been opened, modified or repaired by an unauthorized person.
- If product defects were caused by force majeure such as: lightning strike, electric shock in the electrical network, natural disasters, etc.
- If the defects were caused by damage, improper use or improper transportation.

If the failure was caused by an error in the system to which the product is connected

This guarantee does not change the legal consumer rights valid in EU in relation to those prescribed by the manufacturer.

Authorized service: Feniks-Info d.o.o. Slavonska avenija 72, Zagreb, Croatia servis@feniks-info.hr phone 01/6395 281

Date of purchase: _____

Stamp and signature: _____