
Upute za upotrebu

BEEABeauty

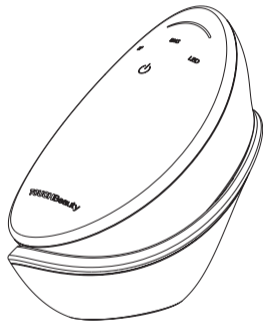TB-1796

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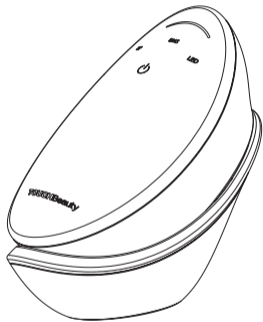
Safety Instructions

TB-1796

To ensure proper operation of this product please read this instruction manual carefully before using the product and keep it for further reference.

WARNING

-  This device is not a toy. Keep out of children's reach.
-  Always verify the device before you start using it. Do not use the device if it is not working properly, appears damaged or broken.
-  Be sure to perform a patch test 24h before using the device for the first time. Apply a small amount of gel to the inner forearm to ensure it does not cause any skin reaction.
-  Avoid using the device for too long on areas you don't want to reduce volume.
-  If used in a dry environment, handling this device may generate static electricity. This is a normal reaction.
-  It does not affect the use of the device.
-  Keep the USB cable and charging dock entrance dry. Do not use the device while it's charging.
-  It is recommended to use the official hydrating gel with this device. Using this device with another gel or cream that contains different ingredients may affect the device or even damage its surface coating.
-  If you have any concerns regarding potential side effects, please consult your doctor before using this device.

 **CAUTION**

Do not use this device in the following situations:

-  Pregnant or nursing women, children; people who have a pacemaker or other electronic medical equipment; if you have skin allergies; people with allergic reactions.
-  If you have active skin diseases (such as wounds, ringworm, eczema, rash, burns, or tissue that is not healing properly, etc.)
-  If you suffer from diabetes, lupus erythematosus, porphyria, congestive heart disease, bleeding disorders, or have an autoimmune disorder (such as HIV or AIDS).
-  If you have an active implant, such as a pacemaker or internal defibrillator; if you suffer from arrhythmia or any other heart condition.
-  If you have any permanent metal implants or piercings in the area of use, or if you have any kind of suspicious lesion in your skin.
-  If you are undergoing radiotherapy or chemotherapy treatments.
-  If you have heat induced diseases such as recurrent herpes.
-  This device is suitable to use on stretch marks, including obesity and pregnancy stretch marks. However, this device does not treat scars.

Do not use this device if you are taking any of the following medications:

- ⊗ If your skin is currently being treated with or has recently been treated with Alpha-Hydroxy Acids (AHAs), Beta-Hydroxy Acids (BHAs), topical isotretinoin, and azelaic acid.
- ⊗ If you have taken any form of Accutane A or Roaccutane A within the past six months.
- ⊗ If you are taking photosensitizing agents or any other doctor-prescribed medication, please check the patient information leaflet before you start using this device. If the patient information leaflet states that it can cause photo-allergic or photo-toxic reactions, or if it specifies to avoid exposure to sunlight, please do not use this device.
- ⊗ If you take anticoagulant medication, including large doses of aspirin, only use the device one week after you stop the treatment.
- ⊗ Use this device on the arms, waist, abdomen, thighs, and buttocks. Never use this device on the face, neck, chest, genital area, and on areas where you have had plastic surgery in the last 3 months.
- ⊗ Do not use light therapy mode if your skin is sensitive to light, prone to rashes or allergic reactions.
- ⊗ Do not use the light therapy if the area of use has infections, eczema, burns, hair follicle inflammation, open wounds, bruises, herpes simplex, lesions, or hematomas.
- ⊗ Do not point or look directly at the light therapy source.
- ⊗ Do not attempt to disassemble, repair, or modify the device by yourself without authorization.
- ⊗ Do not use this device for any other purpose than stated in this instruction manual.

Device Specifications

TB-1796

- Name: REVIVE Multi-Function Body Device BEEA Beauty
- Model No.: TB-1796
- Input: 5V~/0.8A
- Power: 7.4V—/8W
- Type of Battery: Li-ion 7.4V 1000 mAh

● Device Components

Picture 1 Device


Charging Dock

USB cable

Cleaning Cloth

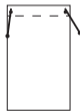Storage Bag


Instrucyion Manual

Device Functions


TB-1760

- ① Power Touch Button
- ② RF Function Touch Button
- ③ EMS Function Touch Button
- ④ LED Function Touch Button
- ⑤ Intensity Level Touch Button
- ⑥ RF/EMS Multipoint Electrodes x 16
- ⑦ RF/EMS Electrode Ring
- ⑧ NIR Lamps x 9
- ⑨ Red Light Lamps x 22
- ⑩ Type-C Charging Terminal

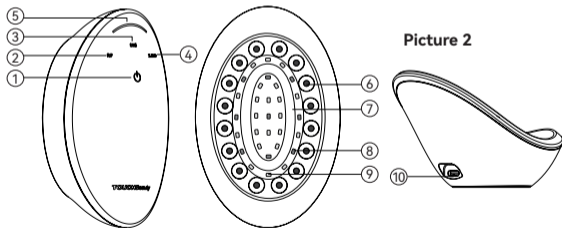※ Charging Instructions

- ☹ This device has a wireless charging system. Place the device on the charging dock. Connect the USB cable to an adapter (adapter not included) and plug the type-c connector into the wireless charging dock terminal. Once the device intensity indicator lights up in a cycle, the device is successfully charging. (Picture 3)
- ☹ Charge this device according to its identification (adapter output: 5V~2.0A). The device does not work while charging. Charging the device takes approximately 2.5 hours. A single charge allows for approximately 60 minutes of use if the battery is full.
- ☹ When the battery is low, the activated function button indicator and the intensity level indicator flash 5 times, and the device shuts down automatically.

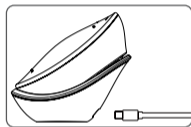Picture 3

※ Turn ON-OFF

- Double-tap the power touch button to turn on the device. By default no function is activated. At this moment, the three function indicators RF/EMS/LED will light up at medium brightness.
- The intensity level indicator does not light up.
- If there is no operation during the standby mode, the device automatically turns off after 1 minute. Double-tap the power touch button at any time to turn on-off the device.

● Mode Switch & Intensity Level Adjustment

- **RF Function Touch Button:** Double-top to turn on-off the RF function.
- **EMS Function Touch Button:** Double-top to turn on-off the EMS function.
- **LED Function Touch Button:** Double-top to turn on-off the LED function.
- **Intensity Level Touch Button:** The intensity level touch button is activated when RF or EMS mode is on. Slide to the right with your finger to increase one intensity level at a time. Slide to the left to decrease one intensity level at a time. You can control the intensity level up to six levels.

How to Use

TB-1796

To ensure the correct operation of this device, please read the instructions carefully before use. Use the device on clean skin. Be sure to remove excess oil and impurities before use.

Treatments & Mode Selection


Treatments	RF	Ring-shaped EMS Microcurrent	LED Red/ Infrared Light	Skincare Product
Stretch Marks/ Wrinkles	●	●	●	Gel
Slimming/Firming	●	●		Body Cream/Gel
Brightening/Rejuvenation	●		●	Gel
Reduce Edema		●		Gel
Repair Sunburns			●	Gel

- Apply an even amount of skincare product to the treatment area before you start using the device. Ensure you apply enough pressure during use to ensure the electrode treatment surface is in full contact with the skin.
- Activate the function you want to use according to the treatment recommendation. Start at a low intensity level and adjust to a level that you feel comfortable with.
- For slimming treatment, it is recommended that you adopt a healthy lifestyle, exercise frequently and give your body enough rest for better results. You can start by applying body cream (optional) to the treatment area. Afterward, apply an even amount of gel to hydrate the skin.
- If you need to purchase the body cream and gel, contact an authorized distributor and purchase these items separately.

Using Techniques

Stretch Marks/Wrinkles

- Use sliding massage movements following the direction of the skin texture for 3 minutes. (Picture 4)
- Use horizontal sliding massage movements for 3 minutes. (Picture 5)
- Use clockwise circular massage movements for 3 minutes. (Picture 6)
- Repeat the massage in counterclockwise movements for 3 minutes. (Picture 7)

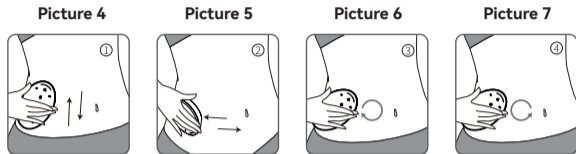
Slimming/Firming

- Waist: Slide the device downwards following the sides of the waist for 10 minutes. (Picture 8)
- Legs: Slide the device upwards along the thigh/calf for 10 minutes. (Picture 9)
- Arms: Slide the device along the upper arm for 10 minutes. (Picture 10)
- Abdomen: Slide the device up and down with a zigzag technique (avoid the navel area) for 10 minutes. (Picture 11)

Picture 8

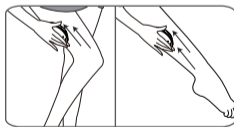Picture 9

Picture 10

Picture 11

- If the skincare product is absorbed by the skin during treatment, reapply it to facilitate the massage movement and improve treatment effectiveness. When you change to a different area of the body, reapply the skincare product to the skin before you start using it.
- The RF function heats up the skin helping to boost collagen regeneration and metabolize subcutaneous fat. If the temperature is not enough, reduce the movement speed; if the temperature is too high, increase the movement speed and adjust it to a comfortable level.
- This device features a motion protection system. During use, make sure you keep moving the device across your skin. If you leave the device in one place without any movement for a long time, the device will stop emitting radiofrequency. When you start moving the device over your skin again, it will start emitting radiofrequency again. However, avoid repetitive movements in the same area for a long time and don't keep the device still in a specific area. Different people have different skin tolerance and sensitivity. Keeping the device without movement on a specific area may cause minor burns. During use, allow the treatment area to have an interval of 2-3 seconds.
- Skin sensation and results may vary from person to person and area of use.
- After the treatment, allow the gel to be completely absorbed by the skin.
- Afterward, you can cleanse the skin.

Using Frequency and duration

- You can choose to activate one or more modes depending on your skin needs. When applying the treatment to the same area, allow at least 1 day rest between treatments.
- It is recommended to use this device 3 times a week. Use the device for a minimum of 6-8 weeks, after this time, use it according to your skin needs. Use the device consistently for better results.

Using duration

- This device has a smart auto timer to give precise treatment cycles. The power will automatically turn off once the treatment time ends. If you want to use the device again, turn on the device and select the function you want to use.

Only 1 function is on RF or EMS or LED: the device turns off automatically after 8 minutes.

Mixed functions RF + EMS or RF + LED or EMS + LED: the device turns off automatically after 10 minutes.

All functions on RF + EMS + LED: the device turns off automatically after 12 minutes.

Cleaning & Maintenance

TB-1796

- 🕒 Regularly clean the device and the accessories to extend its lifetime.
- 🔌 Make sure the device is turned off before cleaning.
- 💧 Always clean the device after use. The waterproof level of this device is IPX6. You can rinse the device with water. However, do not immerse it in water or any other liquid. Wipe off any water or residue with a dry cloth.
- 🚫 Do not use steel wool, abrasive cleansers, chemicals, or corrosive liquids to clean the device.

Storage This device has a lithium battery. Do not place or store the device for long periods in a hot or humid environment. If you do not intend to use the device for a long time, store it after it is fully charged.

Recycling This device has a built-in battery with substances that may pollute the environment. Do not dispose of this device with the household waste and do not throw it into a fire. Please, take this device to a local recycling point for proper disposal. Proper disposal and recycling practices contribute to environmental protection.

Troubleshooting

TB-1796

Problem	Possible Cause	Solution
The device does not work properly.	The battery is running low or out of power.	Charge the device according to the instructions.
You do not feel warm or microcurrent sensation.	Skin is too dry, you did not apply enough gel.	Apply conductive gel.
	You are moving the device too quickly.	Move the device in slower motions.
	The device stopped working.	Keep moving the device on your skin slowly with constant speed.
Note: If the device still does not work properly after using the above solutions, please stop using the device and contact the store or local distributor.		

Warranty

- ※ The warranty period starts from the purchase date.
- ※ The warranty period is 1 year.
- ※ In case of failure, you can return device directly in the store where it was purchased with the purchase invoice, this warranty card, and a detailed description of the defect.

The warranty will be void under the following conditions:

- If the customer does not present the correct warranty card and purchase receipt.
- If the customer did not follow the instructions for using the product
- If the product has been opened, modified or repaired by an unauthorized person.
- If product defects were caused by force majeure such as: lightning strike, electric shock in the electrical network, natural disasters, etc.
- If the defects were caused by damage, improper use or improper transportation.
- If the failure was caused by an error in the system to which the product is connected

This guarantee does not change the legal consumer rights valid in EU in relation to those prescribed by the manufacturer.

Authorized service:

Feniks-Info d.o.o.

Slavonska avenija 72, Zagreb,

servis@feniks-info.hr; tel 01/6395 281

Date of purchase

Stamp and signature