
Upute za upotrebu

BEEABeauty

TB-1959

RF, Sonic & LED
ČEŠALJ ZA NJEGU KOSE

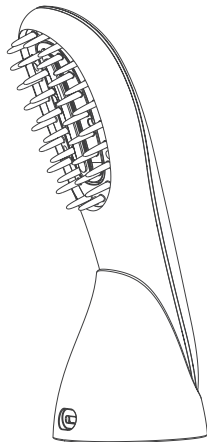Sigurnosne upute

TB-1959

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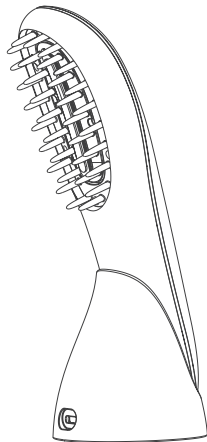
Safety Instructions

TB-1959

To ensure the proper operation of this device, please read this instruction manual carefully before using the device and keep it for further reference.

WARNING

-  Do not wash, rinse or immerse this device in water or any other liquid.
-  This device is not a toy and is not suitable for children. Keep it out of the reach of children.
-  Do not use the device for any purpose other than described in this instruction manual.
-  Keep the comb clean and free of hair styling products, such as mousse, gel, hairspray, etc.
-  If you use this device on dyed hair, it may transfer the color to the comb.
-  This device should only be sent to a service center authorized by the company for maintenance.
-  If handled by unqualified personnel, it may lead to dangerous situations for the user.
-  Use this device carefully by following the specified procedures. Not following the instructions may cause adverse reactions.

CAUTION

-  Do not use this device if you have: heart-related diseases or conditions, facial nerve disorders, diabetes, bleeding disorders, hemophilia, or high blood pressure; if you are undergoing dental treatment; if you have implanted electronic devices; if you have post-pregnancy allergies; if you have implanted metal, plastic, or silicone in your face due to fractures, plastic surgery, etc.; toddlers, children and people that cannot clearly express their wishes; people that have eye sequels, etc.; people that

had eye surgery (if you had a LASIC surgery, please consult your physician before using this device).

- ⊗ Failure to comply may cause hazards, physical discomfort and even blindness.
- ⊗ If you are light-sensitive, consult your physician about the use of this device.
- ⊗ Do not use this device if you have a heat intolerance or your skin does not feel hot temperatures.
- ⊗ Failure to comply may cause burns, electric shock, or injuries.
- ⊗ Do not use this device if: you have allergies, skin diseases/conditions, allergic contact dermatitis, or sensitive skin.
- ⊗ Failure to comply may cause skin and body abnormalities.
- ⊗ Do not use this device if: you are pregnant, in your menstruation cycle, or post-partum.
- ⊗ Failure to comply may cause skin concerns and hormone imbalances.
- ⊗ Do not use this device in the following conditions or areas other than the head: if you feel unwell; while drinking alcohol or taking medicine; if you feel tired; if your skin has wounds, eczema, or is swollen; if you have skin diseases and conditions such as dermatitis, sunburns, etc.; on areas where you have had surgery or recent epilation; on the forehead, lips, eyeballs, upper eyelids, ears, and neck.
- ⊗ Failure to comply may cause skin or physical discomfort and even blindness.
- ⊗ Do not operate this device near flammable or explosive materials, including alcohol, gasoline, thinner, sprays, etc).
- ⊗ Failure to comply may cause a fire and/or electric shock.
- ⊗ Do not throw this device into a fire or heat it. Do not charge, use, or store this device in hot places.
- ⊗ This device has a built-in rechargeable battery that may heat up, catch fire or explode in high temperatures.
- ⊗ If during use your skin becomes hot, red, or swollen caused by the light waves, immediately stop using the device and allow your skin to calm down.
- ⊗ Do not attempt to repair, disassemble, or modify the device to avoid accidents.

Device Specifications


TB-1959

- Name: BEEA Beauty Hair Care Comb RF, Sonic & LED
- Model No.: TB-1959
- Input: 5V /2A
- Power: 3.7V /5W
- Type of Battery: Li-polymer 3.7V 1000mAh

Device Description

TB-1959

- | | |
|-------------------------|--------------------------------------|
| ① Main Unit | ⑦ Laser Mode Button |
| ② Charging Dock | ⑧ Vibration Mode Button |
| ③ Type-C Charging Cable | ⑨ Intensity Level/Charging Indicator |
| ④ Instruction Manual | ⑩ RF Comb Touch-Points |
| ⑤ Power Button | ⑪ LLLT Laser Head |
| ⑥ RF Mode Button | ⑫ Type-C Charging Terminal |

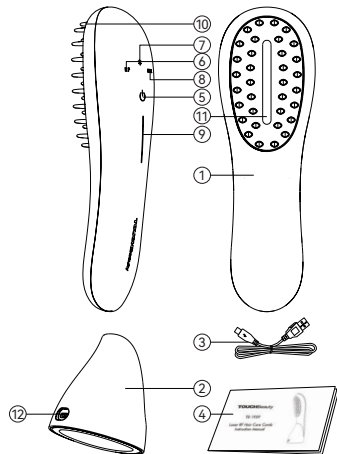
Picture 1

Charging the Device

- Plug the charging cable into the type-c charging terminal of the charging dock and connect the other end of the charging cable to a power supply. Place the main unit on the charging dock to start charging. (Picture 2)
- Charge this device according to its standard identification (adapter output: 5V 2A).
- Charging takes about 5 hours. A full battery allows for approximately 90 minutes of use.
- When charging the device, the white intensity level indicator flashes. When fully charged, the white charging indicator is stable.
- Low Battery Shutdown: All mode button lights flash 5 times with 5 beeping sounds.

Turn ON-OFF

- To Turn ON: Tap the power button twice to turn it ON and enter the standby mode.
- To Turn OFF: With the device turned ON, tap the power button twice to turn it OFF.

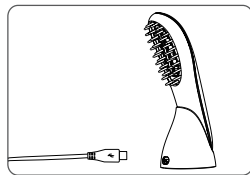Picture 2


How to Use

TB-1959

To ensure the correct operation of this device, please read the instructions carefully before use.

Using the device


- Tap the power button twice to turn the device ON with a sound. Once the device is ON, the RF - laser - vibration mode buttons light up. Select the mode/s according to your needs by tapping twice on the mode you want to use. If you do not select a mode, the device turns OFF automatically after 1 minute. (Picture 3)
- RF Radio Frequency Mode: With the device turned ON, tap on the RF touch button twice to activate the RF mode. The RF mode indicator starts flashing with a beeping sound. Every 8 seconds, the device emits a sound reminding you to move the device to a different area. This mode has 5 intensity levels. You can choose an intensity level that better fits your skin needs and tolerance. (Picture 4)
- Laser Mode: With the device turned ON, tap on the “*” indicator twice to activate the laser mode. The “*” logo starts flashing with a beeping sound. (Picture 5)
- Vibration Massage Mode: With the device turned ON, tap on the SONIC indicator twice to activate the vibration massage mode. The SONIC indicator starts flashing with a beeping sound. Every 4 seconds, the device emits

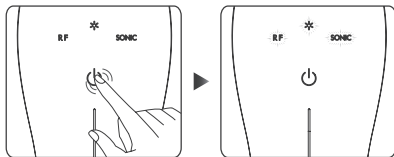
Picture 3

a sound to remind you to change the area of use. This mode has 5 adjustable intensity levels. You can choose an intensity level that better fits your skin needs and tolerance. (Picture 6)


- Intensity Level Adjustment: You can only adjust the intensity level of the RF and vibration mode. Touch and slide the intensity level bar (slide up and down) to adjust the intensity level from levels 1-5. As you slide the intensity bar to adjust the intensity level, the intensity level lights will also increase/decrease with a beeping sound. (Picture 7)
- The device automatically turns OFF after 15 minutes with 3 beeping sounds. It is recommended to use this device once per day. You can use this device with scalp care products to promote absorption.

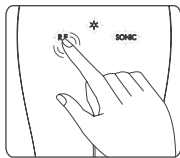
Picture 4

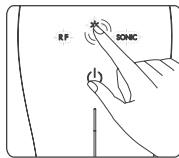Picture 5

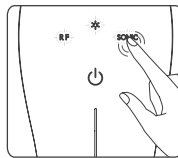Picture 6

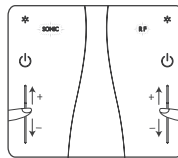Picture 7

Notes

- Do not place this device where it can fall or be pulled into wet areas to prevent external shock.
- Do not clean this device with corrosive liquids to prevent discoloration and damage.
- Store this device out of the reach of children. Operating this device is simple, and the LED light may cause harm or injury to children if they look directly at the light source.
- This device has a lithium battery. Keep the device away from sunlight and store it in a dry location at room temperature (keep it away from hot and humid places). If you don't use the device frequently make sure that you charge it once every 6 months.

Cleaning & Maintenance TB-1959

-  Clean this device with a clean, slightly damp cloth. Do not wash this device with water to prevent damage.
-  Remove any hair or dust from the device and accessories.
-  After cleaning the device with a dampened cloth, place it in a well-ventilated area to dry.

Recycling This device has a built-in battery. If you decide to dispose of this device, please do not throw it in the common household waste. Make sure that you dispose of the battery at a proper recycling center.

Troubleshooting

TB-1959

Problem	Possible Cause	Solution
The device is not working.	The battery is running low or out of power.	Charge the device as guided in this instruction manual.
The device is not charging.	The charging cable is not connected well to the charging dock or power source.	Check the connection between the charging cable and charging terminals to ensure they are well connected.
Note: If the device still does not work properly after using the above solutions, please stop using the device and contact the store or local distributor.		

🔧 Warranty

- ⌘ The warranty period starts from the purchase date.
- ⌘ The warranty period is 1 year.
- ⌘ In case of failure, you can return device directly in the store where it was purchased with the purchase invoice, this warranty card, and a detailed description of the defect.

The warranty will be void under the following conditions:

- If the customer does not present the correct warranty card and purchase receipt.
- If the customer did not follow the instructions for using the product
- If the product has been opened, modified or repaired by an unauthorized person.
- If product defects were caused by force majeure such as: lightning strike, electric shock in the electrical network, natural disasters, etc.
- If the defects were caused by damage, improper use or improper transportation.
- If the failure was caused by an error in the system to which the product is connected

This guarantee does not change the legal consumer rights valid in EU in relation to those prescribed by the manufacturer.

Authorized service:

Feniks-Info d.o.o.

Slavonska avenija 72, Zagreb,

servis@feniks-info.hr; tel 01/6395 281

Date of purchase

Stamp and signature