# SHAPE Shiatsu

## HRV

### UPOZORENJA

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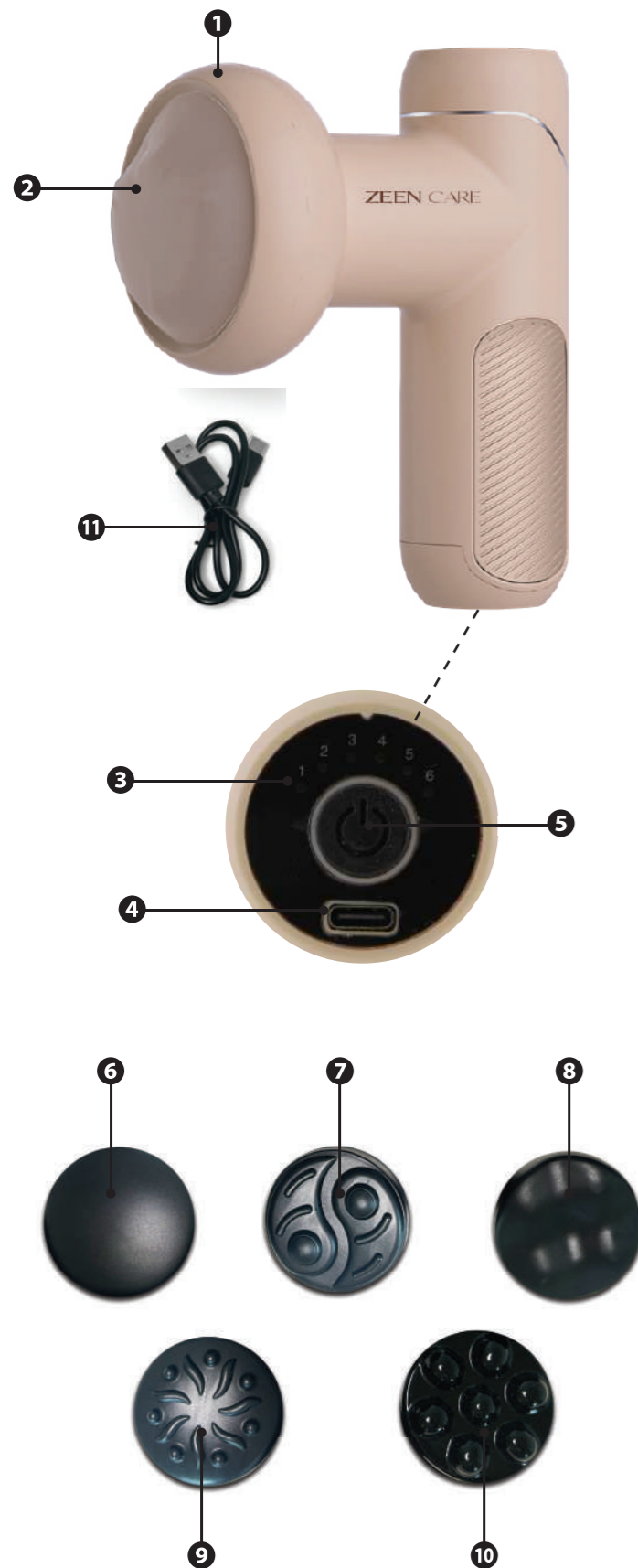
### EN

#### WARNING

- Read all instructions thoroughly and retain for future reference.
- This appliance can be used by children aged 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children should not play with the appliance. Cleaning and maintenance shall not be made by children without supervision.
- Do not immerse product in water.
- For external use only.
- This device must only be used for its intended purpose as described in this manual.
- Do not use on areas that have had cosmetic surgery.
- If you suffer from a skin condition, consult your dermatologist before use.
- This appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- Do not use on cuts or broken skin.
- Keep out of reach from children and pets.
- Do NOT use this appliance if it is not working properly, or if it has been dropped or damaged.
- When charging use a safe temperature range of 0 to 40°C should be used.
- Do not charge with temperature
- Do not charge for longer than 24 hours.
- For hygiene reasons we recommend for single person use only.
- Do not use this appliance on loose skinned parts of your body or any persons showing signs of varicose veins, skin diseases (for example eczema), psoriasis, dermatitis, cracked and abrasive skin, irritation, rashes, male skin tags, skin rash, hypersensitivity or sunburn, blisters, recent scar tissue, other skin conditions. We recommend you consult your physician before treatment if you suffer from any of the above conditions.
- If the power lead is damaged, it must be replaced by the manufacturer or its service agent or similarly qualified persons in order to avoid a hazard.
- Do not use on areas that have had cosmetic surgery.
- If you experience any redness or discomfort stop using immediately.

#### KNOW YOUR UNIT

- |                                 |                       |
|---------------------------------|-----------------------|
| 1.Massager                      | 7.Dual-Point Head     |
| 2.Cover                         | 8.Curved Head         |
| 3.Intensity / Battery Indicator | 9.Circulation Head    |
| 4.USB Charging Socket           | 10.Multi-Point Head   |
| 5.On/Off/ Intensity             | 11.USB Charging Cable |
| 6.Gentle Massage Head           |                       |

#### CHARGING

- 1.Before first use or after a long period of non-use, fully charge the device.
- 2.Connect the USB-C connector of the supplied USB cable to the USB-C port on the device.
- 3.Connect the USB-A end of the cable to a 5 V, min. 1 A USB power adapter (power adapter not included).
- 4.During charging, the intensity/battery indicator may cycle through levels 1 to 6 to display the charging status and battery level.
- 5.Once fully charged, disconnect the device from the power adapter and unplug the cable.

Battery charging time: 2 – 3 hours  
Operating time on full charge: 1 – 2 hours

#### GETTING STARTED

- Press and hold the On/Off/Intensity button for 2 seconds to turn on. The Massager will start at intensity Level 1.
- To increase the intensity, briefly press the On/Off/Intensity button. There are 6 intensity levels.
- To turn off, press and hold the On/Off/Intensity button for 2 seconds.
- Note: The Massager will run for 15 minutes and then automatically switch off.
- Do not apply excessive pressure - simply glide the head over the area. Move slowly using circular motions.
- Note: To minimise soreness it is better to build up to higher intensity levels slowly as your muscles warm up.
- Important: Do not use directly over bones or joints.

#### Total Massage Experience

The SHAPE Shiatsu Massager delivers a powerful, Shiatsu-inspired deep massage designed to support muscle relaxation and overall body comfort. With 6 intensity levels and interchangeable massage heads, you can personalise each session for targeted use across key zones such as the back, shoulders, legs and glutes. USB charging and a compact design enable easy use and effective, customisable massage as part of your daily routine.

#### FITTING THE MESSAGE HEADS

- Ensure the unit is powered off.
- Rotate the Cover anti-clockwise to remove.
- Fit a head in a cavity (5 total), pull firmly to remove.
- Choose your desired massage head, push firmly onto the Massager and reattach the Cover.

#### MESSAGE HEADS

##### Gentle Massage Head

This smooth head is ideal for sensitive areas, delivering the minimum intensity for a very gentle massage treatment.

##### Dual-Point Head

Designed with dual contact points for targeted relief. Ideal for easing muscle knots, trigger points, and deep tension. Best for: Shoulders, back, thighs, and glutes.

##### Curved Head

Features a smooth, flat surface for gentle, relaxing massage. Ideal for sensitive areas or use on sore muscles. Best for: Arms, thighs, glutes, and other large muscle areas.

##### Multi-Point Head

Designed with seven rounded pressure nodes to deliver deep, even stimulation across large muscle groups. Helps relieve muscle tension, improve circulation, and enhance recovery after physical activity.

Best for: Back, thighs, glutes, and other large muscle areas.

##### Circulation Head

Designed with multiple rounded nodes and wave-like ridges to stimulate circulation and support lymphatic flow, providing gentle, full-body massage.

Best for: Arms, legs and abdomen.

#### AFTER USE

- The skin may appear red and slightly irritated after use. This should subside quickly.
- Your muscles may also feel a little sore. This is normal and will lessen with regular use.
- Tips: Ensure you stay hydrated to help flush out any toxins and improve muscle recovery time.

#### CLEANING & MAINTENANCE

1. Use a soft cloth or a gentle cleanser to clean the device. The massage heads can be cleaned directly with water.
2. Keep water, liquids, and humidity away from the device, as moisture may damage the internal electronic components.
3. Store in a cool, dry place. Avoid direct sunlight.

Device Disposal - This device must not be disposed of with household waste. Please take it to a designated collection point for electronic waste in your community or district.

Packaging Disposal - The packaging must not be disposed of with household waste. Please follow the appropriate local procedures for disposal.

#### SPECIFICATIONS

Model: ZC-2456  
Input: 5V == 1A, 5.0W Max  
Battery: Li-Ion 7.4V 1200mAh, 8.88Wh  
Power output: 30 W  
6 speed settings: 1800–3500 RPM  
Battery: 1200 mAh Li-ion  
Runtime: up to 2 hours  
Charging: USB rechargeable  
Rated Voltage: 7.4 V  
Warranty: 12 months

ZEEN CARE