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Do Not Use This Product Under the Following Conditions

- 1.The massager does not operate properly or the massage device is damaged.
- 2.If severe or sudden pain occurs during massage, stop using the product immediately and seek medical assistance.
- 3.If it is confirmed under the above circumstances that the product is the cause of the problem, send the product to a service center for inspection and maintenance

Important Safety Instructions

Before using this product, please read all materials carefully and follow the instructions below:

- 1.Do not place the product in a bathroom or in humid environments.
- 2.Do not leave the massager exposed to direct sunlight or high temperatures for extended periods.
- 3.If the product is used by children, physically or mentally impaired persons, or persons with disabilities, it must be used under adult supervision and assistance.
- 4.This product contains a built-in lithium battery. Please dispose of the product in accordance with applicable environmental protection regulations in your city.

Not Suitable for Use of This Product

If you have any of the following conditions, please consult a physician before use:

- 1.If you are undergoing medical treatment or feel physically unwell
- 2.Patients in the postoperative period
- 3.Persons who consume alcohol
- 4.Minors, adolescents, and children

Use of Rechargeable Battery

- 1.Under normal conditions of use, the service life of the rechargeable battery is approximately 2,500 charge and discharge cycles.
- 2.If the rechargeable battery is not used for a long period, its service life will be shortened. It is recommended to charge the battery at least 1-2 times per month.

Device Disposal

This device must not be disposed of with household waste. Please take it to a designated collection point for electronic waste in your community or district.

Packaging Disposal

The packaging must not be disposed of with household waste. Please follow the appropriate local procedures for disposal.

ZEEN CARE

REVIVO

NECK MASSAGE PILLOW

ZC- 8771

EN

Thank you for choosing this product.

Please read these instructions carefully before use. Keep them in good condition for future reference and information.

Product Effectiveness

- Relieves muscle fatigue and soreness
- Relaxes tension
- Promotes physical and mental relaxation
- Stimulates blood circulation
- Reduces fatigue and exhaustion
- Improves the sub-health condition of the cervical spine

Parameter Description

Model: ZC-8771

Product name: Neck massage pillow

Rated voltage: 3.7 V

Massage duration: approx. 10 minutes

Cover material: velvet, removable and washable

Inner core material: foam

Product dimensions: 25 × 25 × 14 cm

Battery: 1200 mAh (lithium-ion battery)

Power: 3.5 W

Charging voltage: 5 V / 1 A

Function: Forward and reverse rotating massage, red light and heat

Settings: 3 levels

Usage time: approx. 2 hours

Net weight: 495 g

Warranty: 12 months

Charging

- 1.Before first use or after a long period of non-use, fully charge the device.
2. Connect the USB-C connector of the supplied USB cable to the USB-C port on the device.
- 3.Connect the USB-A plug of the cable to a 5 V USB power adapter, minimum 1 A (power adapter not included).
- 4.During charging, the indicator light glows red, and after charging is complete, it turns green.
- 5.After fully charging, disconnect the charger and cable.
- 6.If the charging voltage is unstable, do not use the device while charging. Average charging time is approximately 4 hours. Operating time on a full charge: approximately 2 hours.


Common problems and handling

Serial number	Phenomenon	Possible reason	Approach
1	Massager can't start	Exhausted during standby	It can be started normally after charging.
2	The massager stops working	The timing massage program has been completed or the battery power has been used up	Re-open the battery switch

If there is still an abnormality in the above processing, please contact your local dealer for repair service.




Product configuration

Your U-shaped pillow includes: main unit, charging cable, instruction manual, certificate

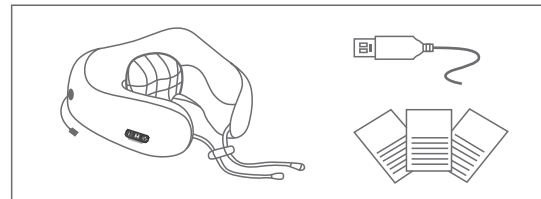
Kneading type

Rope buckle design

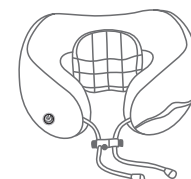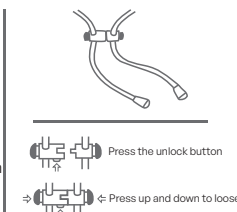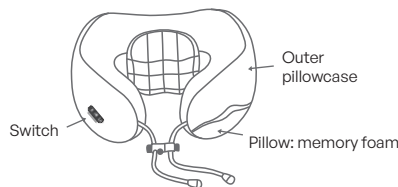One button start

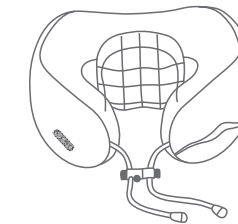One button start

Positive / reverse


Strong bond and weak bond

Red light + heat

Key function operation explanation

Right to left

1. Press and hold the A button for 2 seconds to turn the device ON. Press and hold the A button again for 2 seconds to turn the device OFF. When the device is powered on, it operates in massage mode and is equipped with an automatic shut-off function for safe use.
2. In the power-on working state, press the B button to change the kneading rotation direction. Each press switches the massage direction between Forward and Reverse, allowing a more natural and adjustable massage experience.
3. In the power-on working state, press the C button to select the operating mode. Press once to activate Vibration + Massage. Press again to switch to Massage. Press again to switch to Vibration. This function allows you to choose the preferred mode according to your comfort and needs.
4. In the power-on working state, press the D button to activate the Red light function and heat.
5. The charging indicator features a hidden design. Connect the USB charging cable. After connecting to the power supply, the red indicator light illuminates during charging; when the battery is fully charged, it automatically switches to a green light. The charging time is approximately 240 minutes.

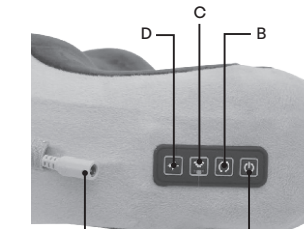USB charging port

A